

PRESSALIT.

Product	PLUS wall track, 15.7", left-grooved, for vertical mounting	
Item number	R9873112	
	VVS:615843134	
	NRF:7501072	
	EAN:5708590369280	
	HMI: 124111	
Color	Anthracite gray	
Product description	Below a selection of products for fitting to the wall track.	
	PLUS wall track, 15.7", left-grooved, for vertical mounting. For height adjustable, left-hand operated support arms. The cover plate hides bolts and wires and facilitates cleaning. Materials: Aluminium and polystyrene. Excl. fixing kit.	
Adjustability	Fixed mounting	
Materials	The wall track is primarily made of 10-micrometer anodized aluminium. The cover plate is made of polystyrene and the ends of the track are capped with ABS polymer.	
Ergonomy	The vertical track enables height adjustment of the different products. If combined with the horizontal track, the available space in the bathroom can be utilized optimally. In environments with multiple users, the vertical track makes it possible to adjust the different products to individual user heights.	
Features	The front cover plate conceals mounting bolts and any cables.	
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.	



PRESSALIT.®

Product	PLUS wall track, 15.7", left-grooved, for vertical mounting		
Item number	R9873112		
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: CaviCide The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.		
	DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.		
Declaration of conformity	This product fulfils all relevant requirements in European Regulation (EU) 2017/745 on medical devices in risk class I. This product is considered to be an accessory for a medical device, as it is not itself a medical device, but is intended to be used together with several particular medical devices to specifically enable the medical devices to be used in accordance with their intended purposes.		



PRESSALIT.

Product				
Item	number			


PLUS wall track, 15.7", left-grooved, for vertical mounting


R9873112

Mounting in different wall types

Concrete >B25


V8610 - Mounting kit, 8 lag bolts $\emptyset 8$ x 70, plugs SX10. To be ordered separately.

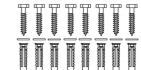


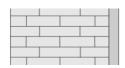

Solid Brick >= Mz12


V8610 - Mounting kit, 8 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.

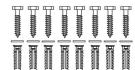

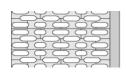

Solid Sand Lime Brick >= Hlz12

V8610 - Mounting kit, 8 lag bolts \emptyset 8 x 70, plugs SX10. To be ordered separately.

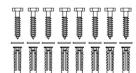





Aerated Concrete P4

V8660 - Mounting kit, 8 countersunk bolts M8 x 35, plugs FTP M 8.

To be ordered separately.

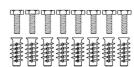


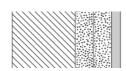

Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood

V8610 - Mounting kit, 8 lag bolts $\emptyset 8$ x 70, plugs SX10. NB! Plugs are not needed.

To be ordered separately.

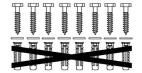



PRESSALIT®

Product Item number	PLUS wall track, 15.7", left-grooved, for vertical mou	nting
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	
	Drywall (2 x 12,5 mm) water-resistant	
	Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.	<u></u>

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com

Pressalit A/S Pressalitvej 1 DK-8680 Ry Tel.: +45 8788 8788 Fax: +45 8788 8789 E-mail: pressalit@pressalit.com