




<p><b>Product</b> <b>Item number</b></p>	<p>PLUS sink bracket with lever control, height adjustable with pneumatic cylinder and manually horizontally adjustable R4680</p>
<p><b>Product description</b></p>	<p>VVS: 654002200 NRF: 6132533 RSK: 7972591 EAN: 5708590309446 HMI: 71757 NBS: N13:30,335 FDA: ILW/D142473</p>  <p>PLUS sink brackets are available in electrically, pneumatically and manually height adjustable versions. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.</p> <p>PLUS sink bracket, pneumatic-assisted. The lever control is placed at the front of the sink. The lever control can also be used as a towel rail and can be placed on either side of the sink. Manually height adjustable by 11.8" (stepless) and adjustable horizontally. Recommendation: Flexible feed and waste system R2084/R2085/R2087. Materials: Aluminium and ABS polymer. Incl. mounting bolts for sink. For mounting on horizontal wall track (not included).</p>
<p><b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b></p>	<p>Vertical and horizontal adjustable 11,82 inch</p> <p>Maximum load: 298 lbs. Suitable for sinks weighing no more than 49 lbs. When the sink is loaded with 298 lbs, each wall fastener will be carrying 66 lbs. <i>Sink brackets need to be able to carry the weight of the sink, as well as support the weight of the person who is leaning against the sink or who is using the sink to rise from seated position.</i> <b>Neither sink nor sink bracket are intended to carry the full body weight of a person, and the test methods applicable to sink brackets reflect this: ISO 17966 prescribes a static test on the front edge of the sink at 75% of maximum user weight, i.e. when the maximum load on a sink bracket is given as 298 lbs and the wash basin weighs 49 lbs, it is – according to ISO 17966 – the equivalent of a user weighing 332 lbs leaning on the front edge of the wash basin.</b></p>



<b>Product</b>	PLUS sink bracket with lever control, height adjustable with pneumatic cylinder and manually horizontally adjustable
<b>Item number</b>	R4680
<b>Materials</b>	The sink bracket is primarily made of 3/8" anodized aluminium. The ends of the bracket cover are capped with ABS polymer. The bracket is available with a lever control consisting of an ABS polymer switch at the end of an anodized aluminium arm
<b>Ergonomy</b>	The height of the sink is easily adjusted the lever control. The sink bracket is compatible with all types of bolt-mounted sink. It's rounded, sleek design makes it easy to clean thoroughly.
<b>Features</b>	The arm that holds the lever control can function as a towel bar. The lever control can be placed on the right- or left-hand side of the sink. This sink bracket is supplied as standard with a pneumatic cylinder for sinks weighing 35.3-48.5 lbs, however, it is possible to order it for sinks weighing 17.6-26.5 lbs or 26.5-35.3 lbs.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

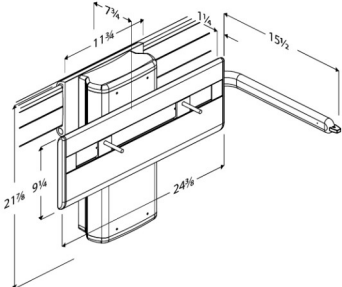

<b>Product</b>	PLUS sink bracket with lever control, height adjustable with pneumatic cylinder and manually horizontally adjustable
<b>Item number</b>	R4680
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



# PRESSALIT®

<b>Product</b>	PLUS sink bracket with lever control, height adjustable with pneumatic cylinder and manually horizontally adjustable
<b>Item number</b>	R4680

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

PRESSALIT®