




PRESSALIT®

Product Item number	PLUS front support, 550 - 620 mm RT874000
Colour Product description	<p>VVS:777539100 NRF:7501042 EAN:5708590314181 HMI: 124195 NBS:N13:426 FDA:INN/D135865</p> <p>White</p>  <p>PLUS front supports are available in three different lengths and colours, they can be used as both front- or back support.</p> <p>PLUS front support, 550 - 620 mm. For PLUS shower seat 450 and PLUS support arms.</p>
Product capacity	Maximum load: 70 kg
Materials Ergonomy Features	<p>Steel and plastics</p> <p>Stabilises and supports users with poor sense of balance.</p> <p>For people with physical disabilities, particularly those with limited upper body strength/stability, or those suffering from involuntary movements or seizures, is the very real risk of falling forwards off, of the toilet. Unlike straps and belts, which can pinch and be uncomfortable against the skin, these adjustable front and back supports attaches directly between two PLUS support arms or between the armrests of a PLUS shower seat. Consequently, users are held securely and comfortably in position, improving their safety and functionality, whilst also assisting carers and offering them greater peace of mind.</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

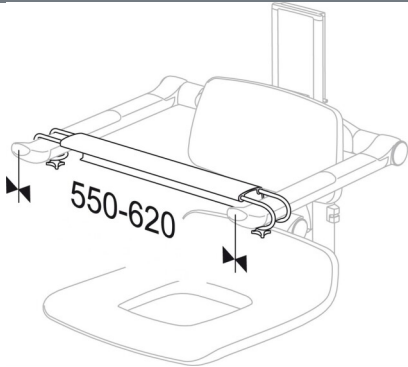

PRESSALIT®

Product	PLUS front support, 550 - 620 mm
Item number	RT874000
Cleaning	<p>PLUS front support and backrest</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>



PRESSALIT®

Product	PLUS front support, 550 - 620 mm
Item number	RT874000

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®