




<b>Product</b> <b>Item number</b>	PLUS front support, 450 - 520 mm RT873000
<b>Colour</b>  <b>Product description</b>	<p>           VVS: 777539000            NRF: 7501039            EAN: 5708590314150            HMI: 124198            NBS: N13:426            FDA: INN/D135865            White         </p>  <p>           PLUS front supports are available in three different lengths and colours, they can be used as both front- or back support.         </p> <p>           PLUS front support, 450 - 520 mm. For PLUS shower seat 310 and PLUS support arms         </p>
<b>Product capacity</b>	Maximum load: 70 kg
<b>Materials</b> <b>Ergonomy</b> <b>Features</b>	<p>Steel and plastics</p> <p>Stabilises and supports users with poor sense of balance.</p> <p>           For people with physical disabilities, particularly those with limited upper body strength/stability, or those suffering from involuntary movements or seizures, is the very real risk of falling forwards off, of the toilet.            Unlike straps and belts, which can pinch and be uncomfortable against the skin, these adjustable front and back supports attaches directly between two PLUS support arms or between the armrests of a PLUS shower seat.            Consequently, users are held securely and comfortably in position, improving their safety and functionality, whilst also assisting carers and offering them greater peace of mind.         </p>
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

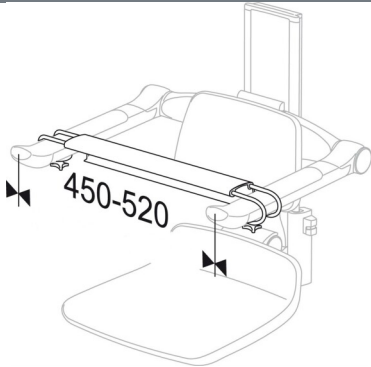

# PRESSALIT®

<b>Product</b>	PLUS front support, 450 - 520 mm
<b>Item number</b>	RT873000
<b>Cleaning</b>	<p>PLUS front support and backrest</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b></p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p><b>NON-DISINFECTANTS:</b></p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b></p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
<b>Declaration of conformity</b>	<p>CE UK MDR This product is in compliance with all relevant requirements of:</p> <p>- Regulation (EU) 2017/745 on medical devices.</p> <p><b>For the UK:</b> - Medical Device Regulation (MDR) 2002 on medical devices.</p>



# PRESSALIT®

<b>Product</b>	PLUS front support, 450 - 520 mm
<b>Item number</b>	RT873000

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®