




<p>Product Item number</p>	<p>PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated R391370112</p>
<p>Color Product description</p>	<p>VVS:615883112 NRF:6191158 RSK:8487453 EAN:5708590368856 HMI: 125893 FDA:ILS/D135866 Anthracite gray</p>  <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions. Suitable for use with toilet pans projecting 21.3/21.7". <i>Replaces discontinued models R3550, RF550, R3565 og RF565.</i></p> <p><i>PLUS fold down grab bar, counter-balanced. Projection 30.3", ergonomically shaped to provide a clear sense of where the fold down grab bar ends. Left hand operated. Height adjustable by 9.5" and adjustable horizontally. Detachable and foldable. Anthracite cover plate for vertical wall track.</i> <i>Materials: The arms are constructed of 10-micrometer anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core</i> <i>For mounting on horizontal wall track (not included).</i></p>
<p>Adjustability Adjustment height Product capacity</p>	<p>Vertical and horizontal adjustable 9,46 inch</p> <p>Maximum load: 330 lbs. When the fold down grab bar is loaded with 330 lbs, each wall fastener in the horizontal wall track is loaded with 97 lbs.</p> <p><i>PLUS fold down grab bars have been tested statically with 1540 lbs (close to the wall), 935 lbs (in the middle) and 495 lbs (the outermost point).</i></p>

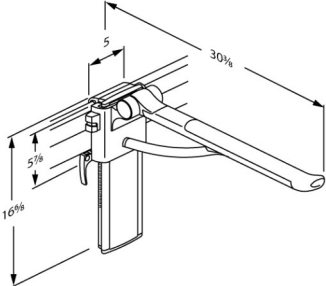

PRESSALIT®

Product	PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated
Item number	R391370112
Materials	<p>The grab bars are primarily constructed in 10-micrometer anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core.</p> <p>The hinge arm is also made of soft SEBS.</p> <p>Internal components and fasteners are primarily made of stainless steel.</p> <p>The hinge barrel is capped with ABS polymer.</p>
Ergonomy	<p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p>
Features	<p>Folds against the wall.</p> <p>The underside is shaped for the mounting of various accessories.</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



Product	PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated
Item number	R391370112
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



Product	PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated
Item number	R391370112
	<p>Intended use The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.</p> <p>Operating environment For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals. For institutional use, the expected life time is 10 years. <i>Life cycle test performed according to ISO 17966:2016.</i></p> <p>For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".</p>