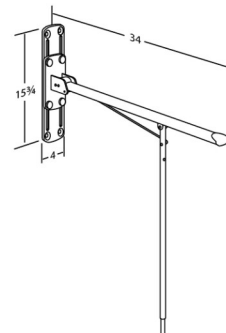




# PRESSALIT®

|                            |   |
|----------------------------|---|
| <b>Product</b>             | VALUE fold down grab bar, height adjustable   |
| <b>Item number</b>         | R1111000  |
| <b>Color</b>               | <p>VVS:615722410<br/>NRF:6191644<br/>RSK:7972413<br/>EAN:5708590176437<br/>HMI: 25414<br/>NBS:N13:436<br/>FDA:ILS/D135866<br/>White</p>   |
| <b>Product description</b> | <p>Fold down grab bar with leg, height adjustable 6".<br/>Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> <p>Fold down grab bar with leg, height adjustable 6".<br/>Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> |
| <b>Adjustability</b>       | Height adjustable   |
| <b>Adjustment height</b>   | 5,91 inch   |
| <b>Product capacity</b>    | <p>Maximum load: 300 lbs.<br/>When the fold down grab bar is loaded with 300 lbs, each mounting bolt is loaded with 69 lbs.</p>   |
| <b>Materials</b>           | <p>Basic material: Steel</p> <p>Surface material: Rilsan</p> <p>End cap: Nylon (PA)</p>   |
| <b>Ergonomy</b>            | <p>Provides firm support when you are sitting, standing, or transferring to and from a toilet.</p> <p>Good, broad support surface for the hand and arm.<br/>The rounded form gives a firm and comfortable grip.</p>   |





# PRESSALIT®

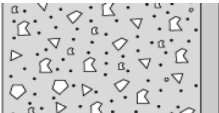
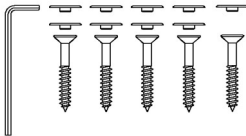

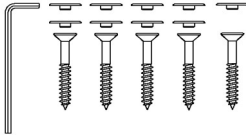
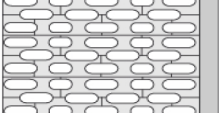
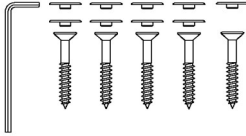
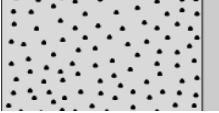
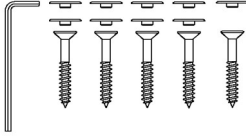
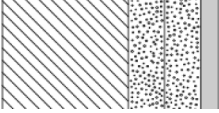
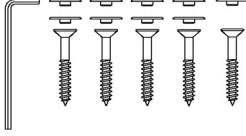
|  |   |
|--|---|
| <b>Product</b><br><b>Item number</b>   | VALUE fold down grab bar, height adjustable<br>R1111000   |
| <b>Features</b>                        | <p>The grab bar can fold up when not in use or to make more space.</p> <p>With leg for extra support.</p> <p>Height adjustable 6".</p>  |
| <b>Warranty</b><br><br><b>Cleaning</b> | <p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p> <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>Clinimax Difficil-S<br/>Erisan OXY+<br/>Sprit (husholdningssprit)<br/>Suma Bac D10<br/>Virkon<br/>Rodalon Desinfektionsspray<br/>Lysoformin Spezial<br/>CaviCide<br/>Sufanios Premium<br/>Tvättsprit 70%<br/>WipeClean UniversalCleaning<br/>Chlorine Disinfection Wipes<br/>Klorilli (Kiilto Clean)</p> <p><b>NON-DISINFECTANTS:</b><br/>Suma Star free D1<br/>Taski 100<br/>Taski 200<br/>Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> |



PRESSALIT®

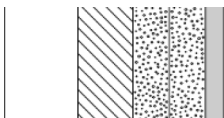
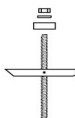


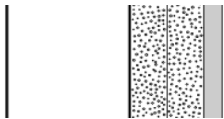
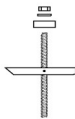

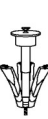
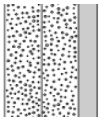

|                    |   |
|--------------------|---|
| <b>Product</b>     | VALUE fold down grab bar, height adjustable |
| <b>Item number</b> | R1111000                                    |

## Mounting in different wall types

|   |   |   |
|---|---|---|
|    | <p><b>Concrete &gt;B25</b></p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included.<br/>Plugs to be ordered separately.</p>   |    |
|   | <p><b>Solid Brick &gt;= Mz12</b></p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included.<br/>Plugs to be ordered separately.</p>                                       |   |
|  | <p><b>Solid Sand Lime Brick &gt;= Hlz12</b></p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included.<br/>Plugs to be ordered separately.</p>                            |  |
|  | <p><b>Aerated Concrete P4</b></p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included.<br/>Plugs to be ordered separately.</p>  |  |
|  | <p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b></p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included.<br/>Plugs to be ordered separately.</p> |  |



# PRESSALIT®

|   |  |   |  |
|---|--|---|--|
| <b>Product</b><br><b>Item number</b>  | VALUE fold down grab bar, height adjustable<br>R1111000  |   |  |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b><br>Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.           |    |  |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b><br>Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar. |    |  |
|  | <b>Drywall (2 x 12,5 mm) water-resistant</b><br>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.   |    |  |

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

[pressalit.com](http://pressalit.com)

Pressalit A/S  
 Pressalitvej 1  
 DK-8680 Ry  
 Tel.: +45 8788 8788  
 Fax: +45 8788 8789  
 E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®