




Product Item number	VALUE support arm, 850 mm, with control panel R1083000
Colour Product description	<p>VVS:615892302 NRF:6191055 RSK:7828177 EAN:5708590346854 HMI: 138195 FDA:ILS/D135866</p> <p>White</p>  <p>VALUE support arm, 850 mm, foldable, incl. operating panel for height adjustment of toilet lifter. Provides firm support when sitting, standing, or transferring to and from a toilet. To be mounted on the right (left hand operated). Height adjustable by 150 mm without using tools.. For fitting to toilet lifter (mounting profile included)</p> <p>VALUE support arm, incl. control panel for TL3/TL4 toilet lifter, 850 mm, foldable, height adjustable 150 mm without the use of tools (left hand operated). Provides firm support when sitting, standing, or transferring to and from a toilet. Materials: Rilsan-coated steel, diecast aluminium and end plug in nylon (PA). For fitting to TL3/TL4 electrical toilet lifter, ALBA wing mounting profile included. For mounting on user's right side.</p>
Adjustability Adjustment height Product capacity	<p>Height adjustable</p> <p>150 mm</p> <p>Maximum load: 120 kg</p>
Materials	<p>The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA). The mounting profile is made from diecast aluminium.</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>



Product	VALUE support arm, 850 mm, with control panel
Item number	R1083000
Cleaning	<p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kilto Clean)</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
WEEE	<p>Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.</p>

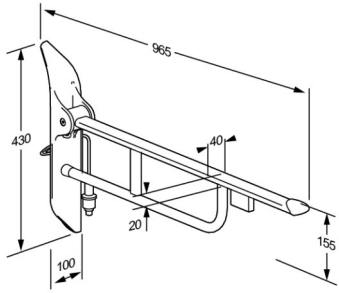

PRESSALIT®

Product	VALUE support arm, 850 mm, with control panel
Item number	R1083000

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®

12. juni 2026