




<p>Product Item number</p>	<p>PLUS fold down grab bar with integrated counter-balance, 33.5", with operating panel, left hand operated R372585108</p>
<p>Color Product description</p>	<p>NRF:6191142 EAN:5708590368672 HMI: 138174 FDA: ILS/D135866</p> <p>Dark blue (RAL 5011)</p>  <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions.</p> <p>PLUS fold down grab bar, with operating panel, counter-balanced, 33.5". Height adjustable 5.9". Foldable. Left handed. Only for use with TL3/TL4 toilet lifter. Mounting profile included. Ergonomically shaped to provide a clear sense of where the support arm ends. Materials: The arms are constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core.</p>
<p>Adjustability Adjustment height Product capacity</p>	<p>Height adjustable 5,91 inch Maximum load: 330 lbs.</p>
<p>Materials Ergonomy</p>	<p>The grab bars are primarily constructed from 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS polymer.</p> <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p>



Product	PLUS fold down grab bar with integrated counter-balance, 33.5", with operating panel, left hand operated
Item number	R372585108
Features	The underside is shaped for the mounting of various accessories.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättspit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>
WEEE	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.

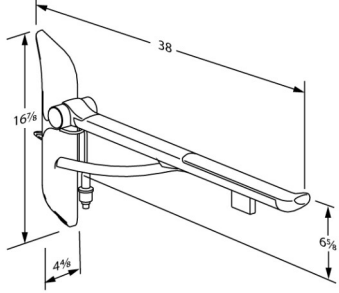

Product	PLUS fold down grab bar with integrated counter-balance, 33.5", with operating panel, left hand operated
Item number	R372585108

Intended use

The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.

Operating environment

For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals.

For institutional use, the expected life time is 10 years.

Life cycle test performed according to ISO 17966:2016.

For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".