




# PRESSALIT®

<b>Product</b>	Shower mattress for MSCT 1 and MSCT 4
<b>Item number</b>	R844721302
	VVS: 777542822 EAN: 5708590373928 HMI: 128434
<b>Colour</b>	Pewter grey
<b>Product description</b>	<p>Shower mattress for MSCT 1 and MSCT 4 shower change trolley, 660 x 2060 mm. Incl. 1.6 m hose.</p> <p>Shower mattress for MSCT 1 and MSCT 4 shower change trolley, 660 x 2060 mm. Incl. 1.6 m hose.</p>
<b>Materials</b>	The mattress cover is welded polyester fabric with PVC-coating. It is UV-resistant and 100% impermeable to water. It contains neither CFC nor latex. The foam filling is flexible polyether-based polyurethane foam, approved in accordance with Oeko-tex 100, Class I.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date.



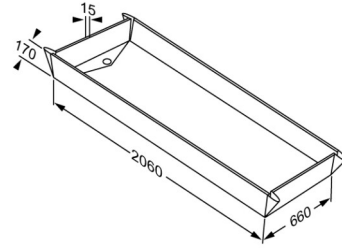

<b>Product</b>	Shower mattress for MSCT 1 and MSCT 4
<b>Item number</b>	R844721302
<b>Cleaning</b>	<p>MSCT mattresses The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Sprit (husholdningssprit) Tvättspirit 70% Safe Zone Plus Selgiene Extreme Microsan Accel PREvention Wipes Lingettes</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	<p>MDR-ACC This product fulfils all relevant requirements in Regulation (EU) 2017/745 on medical devices in risk class I. This product is considered to be an <b>accessory</b> for a medical device, as it is not itself a medical device, but is intended to be used together with several particular medical devices to specifically enable the medical devices to be used in accordance with their intended purposes.</p>