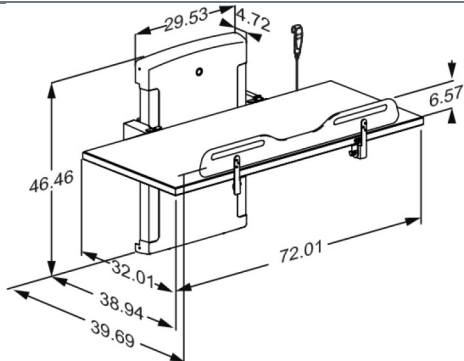




PRESSALIT®

<b>Product</b> <b>Item number</b>	CT 4100 change table, electrically foldable, electrically height adjustable R8595572000
<b>Color</b>	<div> <div> EAN:5708590356952  FDA:INQ/D124319 </div> <div>  </div> </div> White
<b>Product description</b>	<p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>The CT 4000 series has been designed for public changing facilities and schools.</p> <p>IAPMO-certified.  CT 4100 change table, electrically foldable.  Electrically height adjustable 27.6"; from 11.8" - 39.4", with buttons placed in table front.  Incl. safety rail, body belt and wired remote control.  Lying surface 32" x 72".  Outer measurements 39.7" x 72".  Excl. mounting kit.</p>
<b>Adjustability</b>	Height adjustable
<b>Adjustment height</b>	27,58 inch
<b>Product capacity</b>	<p>Maximum user weight (MUM): 500 lbs.  Maximum load (SWL): 511 lbs. (User weight + accessories)  When the table is loaded with 500 lbs, each mounting bolt supports 83 lbs.  <i>This product has been tested statically with a load of 1279 lbs as the IEC 60601-1 standard requires testing with a safety factor 2.5.</i></p>
<b>Speed up</b>	<p>Speed up: 0,59 inch/s  Speed down: 0,79 inch/s</p>
<b>Noise level</b>	48 dB



# PRESSALIT®

<b>Product</b> <b>Item number</b>	CT 4100 change table, electrically foldable, electrically height adjustable R8595572000
<b>Materials</b>  <b>Ergonomy</b>	<p>Powder-laquered stainless steel, lying surface in Compact Density Fiberboards (CDF).</p> <p>The ergonomically designed safety rail has a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user. Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.</p>
<b>Mains connection (V)</b>	120
<b>Actuator (motor)</b>	LA40 Careline (Linak) 4000 N, IPX6
<b>Control box</b>	CA63 Careline (Linak), IPX6
<b>Wired remote control</b>	HB30 (Linak) IP66
<b>Power consumption (maximum load)</b>	Power consumption by maximum load: 420 W Power consumption by standby: 2 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date.



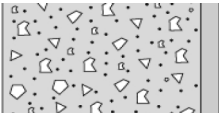
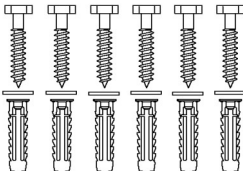

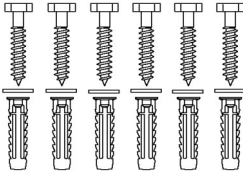
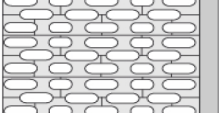
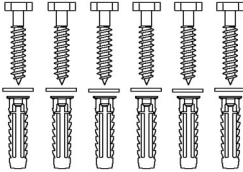
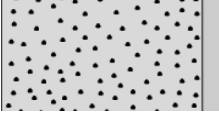
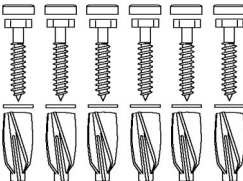
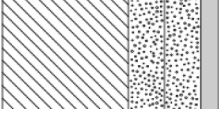
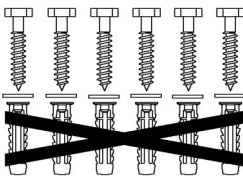
<b>Product</b>	CT 4100 change table, electrically foldable, electrically height adjustable
<b>Item number</b>	R8595572000
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes OxiGenesis from ENVIROX</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>



PRESSALIT®

<b>Product</b>	CT 4100 change table, electrically foldable, electrically height adjustable
<b>Item number</b>	R8595572000

## Mounting in different wall types

	<p><b>Concrete &gt;B25</b></p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Solid Brick &gt;= Mz12</b></p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Solid Sand Lime Brick &gt;= Hlz12</b></p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p><b>Aerated Concrete P4</b></p> <p>V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b></p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	



# PRESSALIT®

<b>Product</b> <b>Item number</b>	CT 4100 change table, electrically foldable, electrically height adjustable R8595572000
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b></p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b></p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p><b>Drywall (2 x 12,5 mm) water-resistant</b></p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

[pressalit.com](http://pressalit.com)

Pressalit A/S  
 Pressalitvej 1  
 DK-8680 Ry  
 Tel.: +45 8788 8788  
 Fax: +45 8788 8789  
 E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®