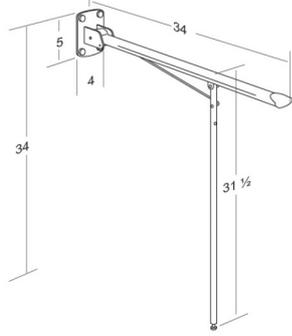




Product Item number	VALUE fold down grab bar, fixed height R1102000
Color Product description	<p>VVS:615722380 NRF:6191677 RSK:7972411 EAN:5708590192345 HMI: 29695 NBS:N13:436 DB: 2109684 FDA:ILS/D135866</p> <p>White</p>  <p>VALUE fold down grab bar with leg, fixed, foldable, length 33.5", height 31.5". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, pressure-cast aluminium, plastic. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p>
Adjustability Product capacity	<p>Fixed mounting</p> <p>Maximum load: 300 lbs. When the fold down grab bar is loaded with 300 lbs, each mounting bolt is loaded with 180 lbs.</p>
Materials Ergonomy	<p>Basic material: Steel</p> <p>Surface material: Rilsan</p> <p>End cap: Nylon (PA)</p> <p>Provides firm support when you are sitting, standing, or transferring to and from a toilet.</p> <p>Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.</p>

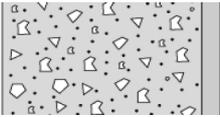
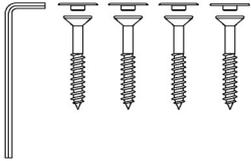
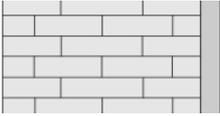
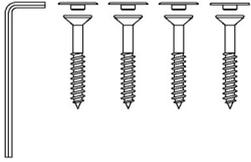
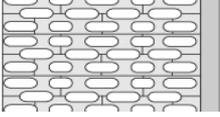
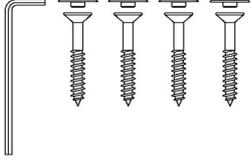
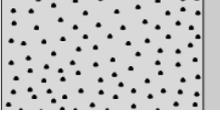
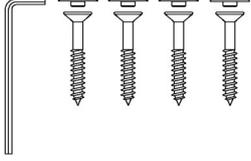
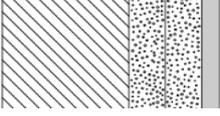
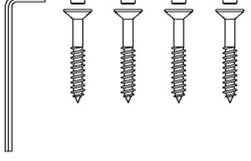


Product	VALUE fold down grab bar, fixed height
Item number	R1102000
Features	<p>The grab bar can fold up when not in use or to make more space.</p> <p>With leg for extra support.</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean)</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>

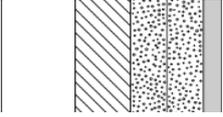
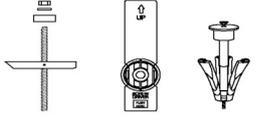
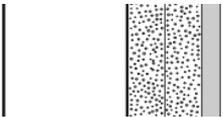
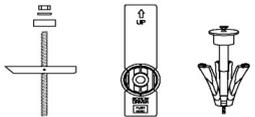
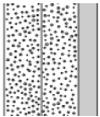


Product	VALUE fold down grab bar, fixed height
Item number	R1102000

Mounting in different wall types

	Concrete >B25 V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Brick >= Mz12 V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Aerated Concrete P4 V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	



Product Item number	VALUE fold down grab bar, fixed height R1102000
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or Wingt anchors or similar. 
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or Wingt anchors or similar. 
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. 

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.