




| | |
|---|--|
| <p>Product Item number</p> | <p>PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated R372370000</p> |
| <p>Color Product description</p> | <p>VVS:615877100 NRF:6191124 EAN:5708590368511 HMI: 138189 FDA:ILS/D135866 White</p>  <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions.</p> <p>PLUS fold down grab bar, counter-balanced, 27.6". Height adjustable 5.9". Foldable. Left handed. Only for use with TL3/TL4 toilet lifter. Mounting profile included. Ergonomically shaped to provide a clear sense of where the support arm ends. Materials: The arms are constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core.</p> |
| <p>Adjustability Adjustment height Product capacity</p> | <p>Height adjustable 5,91 inch Maximum load: 330 lbs.</p> |
| <p>Materials Ergonomy</p> | <p>The grab bars are primarily constructed from 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS polymer.</p> <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p> |



| | |
|--------------------|--|
| Product | PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated |
| Item number | R372370000 |
| Features | The underside is shaped for the mounting of various accessories. |
| Warranty | Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. |
| Cleaning | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |

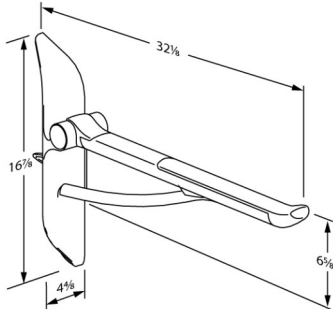

| | |
|--------------------|--|
| Product | PLUS fold down grab bar with integrated counter-balance, 27.6", left hand operated |
| Item number | R372370000 |

Intended use

The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.

Operating environment

For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals.

For institutional use, the expected life time is 10 years.

Life cycle test performed according to ISO 17966:2016.

For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".