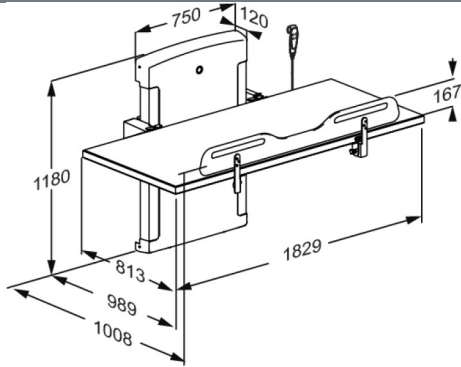




PRESSALIT®

Product Item number	CT 4000 change table, manually foldable, electrically height adjustable R8584572000
Colour	VVS: 777541100 RSK: 8565085 EAN: 5708590358130 HMI: 125834 White
Product description	 Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user. The CT 4000 series has been designed for public changing facilities and schools. CT 4000 change table, manually foldable. Electrically height adjustable 700 mm; from 300 mm -1000 mm, with buttons placed in table front. Incl. safety rail, body belt and wired remote control. Lying surface 813 x 1829 mm in compact laminate. Outer measurements 1008 mm x 1829 mm. Excl. mounting kit.
Adjustability Adjustment height Product capacity Speed upwards Noise level	Height adjustable 700 mm Maximum user weight (MUM): 227 kg. Safe working load (SWL): 232 kg. When the table is loaded with 227 kg, each mounting bolt supports 37,5 kg <i>This product has been tested statically with a load of 580 kg as the IEC 60601-1 standard requires testing with safety factor 2.5.</i> Speed up mm: 15 mm/s Speed down mm: 20 mm/s 48 dB



Product	CT 4000 change table, manually foldable, electrically height adjustable
Item number	R8584572000
Materials	Powder-laquered stainless steel, lying surface in Compact Density Fiberboards (CDF).
Ergonomy	The ergonomically designed safety rail has a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user. Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.
Features	• Built-in switches in the worktop – for easy adjustment of the height. You can stay close to the user who is lying on the worktop. The built-in switches are placed conveniently close to the helper. It is also possible to lower and raise the table by using the wired hand control. - Lifting height up to 1000 mm to provide proper and safe working height for the helper. • The mattress is reversible so it is possible to turn it over, if one side is worn out. The reversible mattress is fixed to the table so it will not move when changing a restless person or while folding/unfolding the table.
Mains connection (V)	230
Actuator (motor)	LA40 Careline (Linak) 4000 N, IPX6
Control box	CA63 Careline (Linak), IPX6
Wired remote control	HB30 (Linak) IP66
Power consumption (maximum load)	Power consumption by maximum load, watt: 420 W Power consumption by standby: 2 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



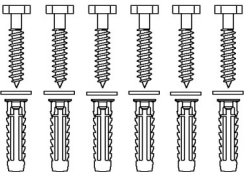
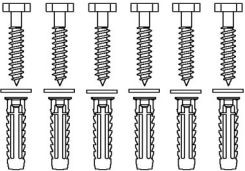
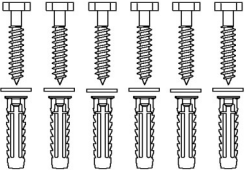
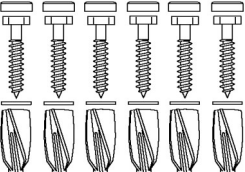
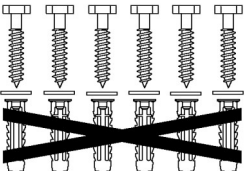
PRESSALIT®

Product	CT 4000 change table, manually foldable, electrically height adjustable
Item number	R8584572000
Cleaning	The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes OxiGenesis from ENVIROX NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.
WEEE	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.



Product	CT 4000 change table, manually foldable, electrically height adjustable
Item number	R8584572000

Installation in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Aerated Concrete P4 V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.	



Product Item number	CT 4000 change table, manually foldable, electrically height adjustable R8584572000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	  	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	  	
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.