




PRESSALIT®

Product Item number	PLUS support arm with integrated counter-balance, 700 mm, with operating panel, left hand operated R372670112
Colour	Anthracite grey
Product description	VVS:615878312 NRF:1617117 EAN:5708590391113  All PLUS support arms come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions. PLUS support arm with operating panel, counter-balanced, 700 mm. Height adjustable 150 mm. Foldable. Left handed. Only for use with TL5 toilet lifter equipped with toilet pan projecting 540/550 mm. Mounting profile included. Ergonomically shaped to provide a clear sense of where the support arm ends. Materials: The arms are constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core.
Adjustability	Height adjustable
Adjustment height	150 mm
Product capacity	Maximum load: 150 kg
Materials	The arms are primarily constructed in 10-micrometre anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS polymer.
Ergonomy	The oval shape gives an easy and comfortable surface to lean and push against. The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends. Soft and slip-resistant material that has a warm feel.

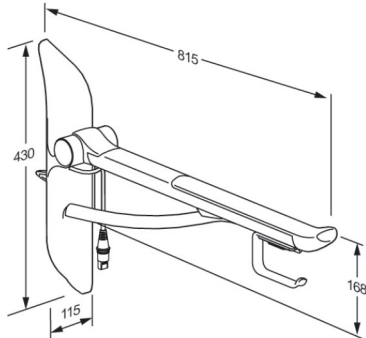

PRESSALIT®

Product	PLUS support arm with integrated counter-balance, 700 mm, with operating panel, left hand operated
Item number	R372670112
Features	The underside is shaped for the mounting of various accessories.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättssprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate. DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.
WEEE	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.



PRESSALIT®

Product	PLUS support arm with integrated counter-balance, 700 mm, with operating panel, left hand operated
Item number	R372670112

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®

8. juli 2026