



Toilet seat Dania without cover, 100 mm raised
With stabilising buffers and special hinges adapted to raised toilet
seat. Incl. splash guard.

Item number

R24000
 VVS:615686000
 NRF:6190527
 EAN:5708590070483
 HMI: 1815
 NBS:N13:300

Toilet seat, Dania, with stabilising buffers, w/o cover,
 raised 100 mm, incl. splashguard. With raised cross bar
 hinge





Colour

White

Description

Toilet seat, Dania, with stabilising buffers, w/o cover,
 raised 100 mm, incl. splashguard
 With raised cross bar hinge

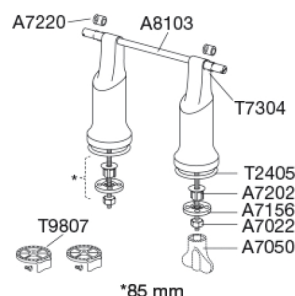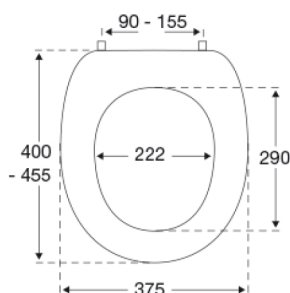

Materials

TOILET SEAT:

The material is colour ingrained duroplast (UF A 10 = urea formaldehyde) that contains no environmentally hazardous substances. UF plastic comprises approx. 70% urea formaldehyde resin, 30% cellulose, as well as 5% minerals, pigments, lubricants and moisture content. Free formaldehyde is released during the production process. Pressalit toilet seats have been measured after curing ($< 0.01 \text{ mg/m}^3$) to emit significantly below the permitted EU limit of 0.08 mg/m^3 .

BUFFERS:

EVA (copolymer made of ethylene and vinyl acetate).

MOUNT/HINGES:

Stainless steel and thermoplastic.

Product capacity

This product has been tested statically up to 225 kg. As the European standard for technical aids for disabled, BS/EN 12182:2012, requires a static test with safety factor 1.5, the nominal maximum load for this product is therefore 150 kg.

Weight

2,8 kg

Gross weight

3,4 kg



PRESSALIT®

Toilet seat Dania without cover, 100 mm raised
With stabilising buffers and special hinges adapted to raised toilet
seat. Incl. splash guard.

Weight per crate 0 kg

Qty per crate 0

Qty per pallet 0

Cleaning Use a mild soap solution to clean the toilet seat. Please take care that the seat and hinges are not left damp. For disinfecting the seat, e.g. ethanol 85% can be used. Make sure that the product is wiped down with a soft rag wrung out in clear water and dried off with a soft, dry cloth immediately after disinfection.

Avoid contact of the seat and hinges with abrasive, corrosive or chlorine-based cleaners, as these may cause damage or lead to flash rust. Therefore, when cleaning the bowl, make sure that the seat and cover are kept in an upright position until all cleaning agent has been flushed away.

Warranty Your Pressalit toilet seat is subject to a 10-year product warranty against presentation of receipt. You should contact the dealer or store where you purchased the seat in the event of a complaint. In the event of complaints about defective materials and/or manufacturing defects and/or functional limitations ascertained within a 10-year period, following proper use and cleaning of the toilet seat (we reserve the right to examine the toilet seat at our factory), Pressalit will deliver spare parts or replace the toilet seat by an identical or similar model. If the color is no longer available, Pressalit reserves its right to replace the toilet seat by a similar model in the color white. Upon replacement or delivery of a spare part, the warranty period will not be extended nor renewed. If replacement is not possible, Pressalit may instead refund the purchase price.

Subject to changes in the product range, technical changes and printing errors.

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®