




<p>Product Item number</p>	<p>MSCT 4 shower change trolley, manually height adjustable R9504717000</p>
<p>Colour Product description</p>	<p>EAN:5708590389844 FDA:INW/D546021</p>  <p>White</p> <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>MSCT 4 shower change trolley, 700 x 1700 mm, four fold-down sides. Manually height adjustable with foot pedal (hydraulic). 500 mm lifting range, from 550 to 1050 mm.</p> <p>Mattress to be ordered separately.</p>
<p>Adjustability Adjustment height Product capacity</p>	<p>Height adjustable</p> <p>500 mm</p> <p>Maximum user weight (MUM): 200 kg. Maximum load (SWL): 250 kg. (weight + water + accessories). <i>This product has been tested statically with a load of 375 kg as the ISO 17966 standard requires testing with safety factor 1.5.</i></p>
<p>Materials</p>	<p>Plastics, compact laminate and powder-lacquered stainless steel.</p>



Product	MSCT 4 shower change trolley, manually height adjustable
Item number	R9504717000
Ergonomy	<ul style="list-style-type: none">• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user.• One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table.• The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training.• Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.
Features	<ul style="list-style-type: none">• The trolley can be moved from room to room for easy transport and multiple users.• Intuitive, easy-to-use drop down rails with one hand operation.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.

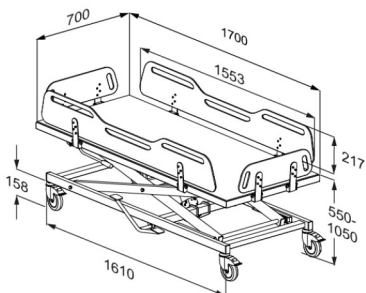

Product	MSCT 4 shower change trolley, manually height adjustable
Item number	R9504717000
Cleaning	<p>MSCT 1, MCT 2, MSCT 4 The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
Declaration of conformity	<p>CE UK MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices. For the UK: - Medical Device Regulation (MDR) 2002 on medical devices.</p>