




<p>Product Item number</p>	<p>VALUE mounting plate, manually height adjustable R1214000</p>
<p>Colour Product description</p>	<p>VVS:654011410 NRF:6191676 EAN:5708590176499 HMI: 25578 NBS:N13:30,335 FDA:ILW/D142473</p> <p>White</p>  <p>Mounting plate incl. brackets for wallhung wash basin, height adjustable 150 mm (recommendation R2065/R2067 (UK: R2075/R2077) flexible feed/waste system). Incl. mounting kit for brick wall</p> <p>Mounting plate incl. brackets for wallhung wash basin, height adjustable 150 mm. (recommendation R2065/R2067 (UK: R2075/R2077) flexible feed/waste system). Incl. mounting kit for brick wall</p>
<p>Adjustability Adjustment height Product capacity</p>	<p>Height adjustable</p> <p>150 mm</p> <p>Maximum load: 80 kg. Suitable for wash basins weighing no more than 22 kg. When the wash basin is loaded with 80 kg, each wall fastener will be carrying 22 kg. <i>Wash basin brackets need to be able to carry the weight of the wash basin, as well as support the weight of the person who is leaning against the wash basin or who is using the wash basin to rise from seated position.</i> Neither wash basin nor wash basin bracket are intended to carry the full body weight of a person, and the test methods applicable to wash basin brackets reflect this: ISO 17966 prescribes a static test on the front edge of the wash basin at 75% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 80 kg and the wash basin weighs 22 kg, it is – according to ISO 17966 – the equivalent of a user weighing 77 kg leaning on the front edge of the wash basin.</p>
<p>Materials</p>	<p>Basic material: Steel</p> <p>Surface material: Rilsan</p>



Product	VALUE mounting plate, manually height adjustable
Item number	R1214000
Features	Height adjustable 150 mm.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>

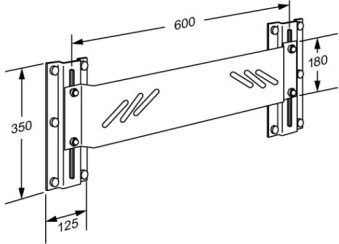

PRESSALIT®

Product	VALUE mounting plate, manually height adjustable
Item number	R1214000

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®