




Product Item number	VALUE support arm, 850 mm R1085000
Colour Product description	<p>VVS:615892502 NRF:6191056 RSK:7828178 EAN:5708590346861 HMI: 138199 FDA:ILS/D135866</p> <p>White</p>  <p>VALUE support arm, 850 mm, for toilet lifter, foldable. Provides firm support when sitting, standing, or transferring to and from a toilet. To be mounted on the right (left hand operated). Height adjustable by 150 mm without using tools. For fitting to toilet lifter (mounting profile included).</p> <p>VALUE support arm, 850 mm, foldable, height adjustable 150 mm without the use of tools (left hand operated). Provides firm support when sitting, standing, or transferring to and from a toilet. Materials: Rilsan-coated steel, diecast aluminium and end plug in nylon (PA). For fitting to TL3/TL4 manual toilet lifter, ALBA wing mounting profile included. For mounting on user's right side.</p>
Adjustability Adjustment height Product capacity	<p>Height adjustable</p> <p>150 mm</p> <p>Maximum load: 120 kg</p>
Materials	<p>The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA). The mounting profile is made from diecast aluminium.</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>

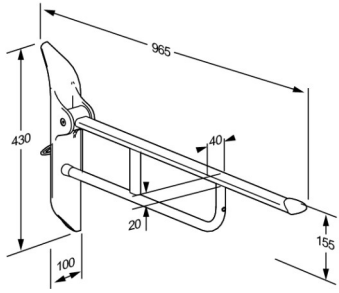

Product	VALUE support arm, 850 mm
Item number	R1085000
Cleaning	<p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kilto Clean)</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>