




PRESSALIT®

Product Item number	VALUE support arm, 850 mm, with operating panel, left hand operated R102685000
Colour	VVS:615892602 NRF:1617150 EAN:5708590391069 White
Product description	 VALUE and VALUE II support arms with mounting profile are specifically designed for use with SELECT TL5 toilet lifters. The height adjustment of the toilet lifter is done only with the control panel on the support arm. VALUE support arm with operating panel, 850 mm. Height adjustable 150 mm. Foldable. Left hand operated. Only for use with TL5 toilet lifter. Mounting profile included. For mounting on user's right side. Materials: Rilsan-coated steel, diecast aluminium and end plug in nylon (PA). Excl. mounting kit.
Adjustability	Height adjustable
Adjustment height	150 mm
Product capacity	Maximum load: 120 kg
Materials	The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA). The mounting profile is made from diecast aluminium.
Ergonomy	Provides firm support when sitting, standing, or transferring to and from a toilet.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

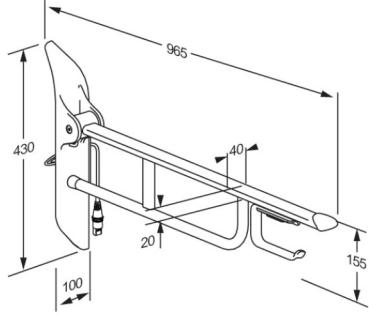

Product	VALUE support arm, 850 mm, with operating panel, left hand operated
Item number	R102685000
Cleaning	The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged. For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean) NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.
WEEE	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.



PRESSALIT®

Product	VALUE support arm, 850 mm, with operating panel, left hand operated
Item number	R102685000

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®

8. juli 2026