




<p><b>Product</b> <b>Item number</b></p>	<p>MCT 3 change trolley, electrically height adjustable R9612717000</p>
<p><b>Colour</b>  <b>Product description</b></p>	<p>EAN:5708590389899</p>  <p>White</p> <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>MCT 3 change trolley, 700 x 1700 mm, two fold-down sides. Electrically height adjustable with wired remote control. 500 mm lifting range, from 450 to 950 mm.</p> <p>Mattress to be ordered separately.</p>
<p><b>Adjustability</b> <b>Adjustment height</b>  <b>Product capacity</b>  <b>Speed upwards</b>  <b>Noise level</b></p>	<p>Height adjustable</p> <p>500 mm</p> <p>Maximum user weight (MUM): 200 kg. Maximum load (SWL): 205 kg. <i>This product has been tested statically with a load of 512,5 kg as the IEC 60601-1 standard requires testing with safety factor 2.5.</i></p> <p>Speed up mm: 17 mm/s Speed down mm: 17 mm/s</p> <p>50 dB</p>
<p><b>Materials</b></p>	<p>Plastics, compact laminate and powder-lacquered steel.</p>



<b>Product</b>	MCT 3 change trolley, electrically height adjustable
<b>Item number</b>	R9612717000
<b>Ergonomy</b>	<ul style="list-style-type: none"><li>• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user.</li><li>• One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table.</li><li>• The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training.</li><li>• Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.</li></ul>
<b>Features</b>	<ul style="list-style-type: none"><li>• The trolley can be moved from room to room for easy transport and multiple users.</li><li>• Intuitive, easy-to-use drop down rails with one hand operation.</li></ul>
<b>Mains connection (V)</b>	230
<b>Actuator (motor)</b>	LA40 (Linak) IPX6
<b>Control box</b>	CA30 (Linak) IPX6
<b>Wired remote control</b>	HB30 (Linak) IPX6
<b>Power consumption (maximum load)</b>	Power consumption by maximum load, watt: 280 W Power consumption by standby: 0,5 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.

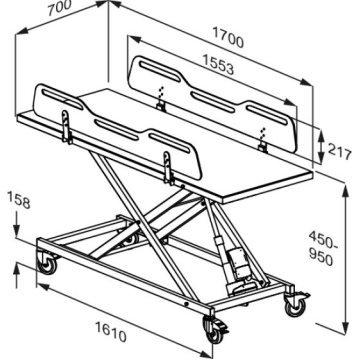

<b>Product</b>	MCT 3 change trolley, electrically height adjustable
<b>Item number</b>	R9612717000
<b>Cleaning</b>	<p>MCT 3, MCT 5</p> <p>The product is to be cleaned with a damp cloth that has been wrung out in a mild soap mixture, such as universal cleaning solution. After cleaning, it is essential to wipe the product with a dry, soft cloth. This step is crucial to prevent the materials from discoloration. Additionally, it helps prevent any potential degradation of the surface paint,</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>WEEE</b>	<p>Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.</p>
<b>Declaration of conformity</b>	<p>CE UK MDR-MASK-RoHS This product is in compliance with all relevant requirements of:</p> <ul style="list-style-type: none"><li>- Regulation (EU) 2017/745 on medical devices</li><li>- Council Directive 2006/42/EC on machinery</li><li>- Council Directive 2011/65/EC on RoHS. <b>For the UK:</b></li><li>- Medical Device Regulation (MDR) 2002 on medical devices.</li><li>- Supply of Machinery (Safety) Regulations 2008</li><li>- Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012.</li></ul>



# PRESSALIT®

<b>Product</b>	MCT 3 change trolley, electrically height adjustable
<b>Item number</b>	R9612717000

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

PRESSALIT®