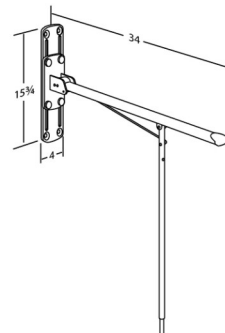




PRESSALIT®

Product	VALUE fold down grab bar, height adjustable
Item number	R1111000
Color	VVS:615722410 NRF:6191644 RSK:7972413 EAN:5708590176437 HMI: 25414 NBS:N13:436 FDA:ILS/D135866 White
Product description	Fold down grab bar with leg, height adjustable 6". Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately. Fold down grab bar with leg, height adjustable 6". Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.
Adjustability	Height adjustable
Adjustment height	5,91 inch
Product capacity	Maximum load: 300 lbs. When the fold down grab bar is loaded with 300 lbs, each mounting bolt is loaded with 69 lbs.
Materials	Basic material: Steel Surface material: Rilsan End cap: Nylon (PA)
Ergonomy	Provides firm support when you are sitting, standing, or transferring to and from a toilet. Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.





PRESSALIT®

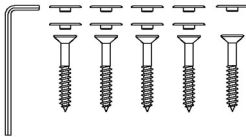
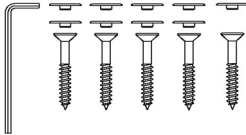
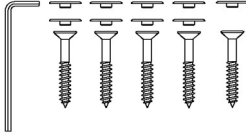
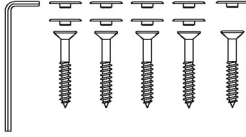
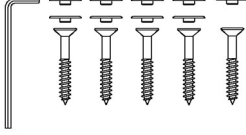
Product Item number	VALUE fold down grab bar, height adjustable R1111000
Features	The grab bar can fold up when not in use or to make more space. With leg for extra support. Height adjustable 6".
Warranty Cleaning	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged. For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean) NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.



PRESSALIT®

Product	VALUE fold down grab bar, height adjustable
Item number	R1111000

Mounting in different wall types

	Concrete >B25 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Brick >= Mz12 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Aerated Concrete P4 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	



PRESSALIT®

Product Item number	VALUE fold down grab bar, height adjustable R1111000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®