
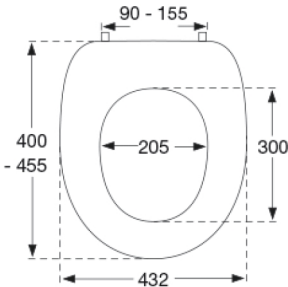
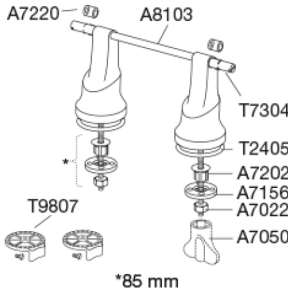
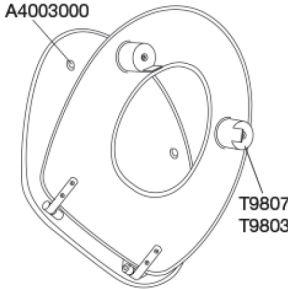
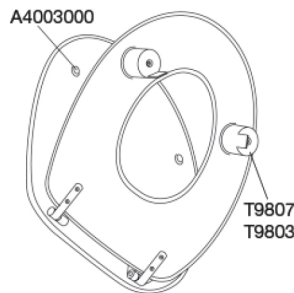
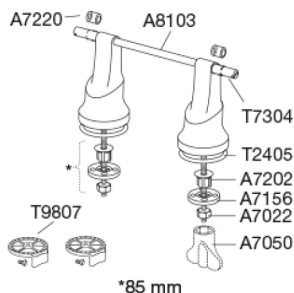
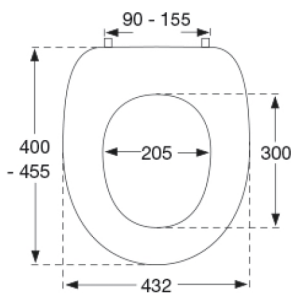




# PRESSALIT®

**Toilet seat Ergosit with cover, 50 mm raised**  
**With stabilising buffers and special hinges adapted to raised toilet seat.**

<b>Item number</b>	R20000 VVS:615691000 NRF:6190519 EAN:5708590094779 HMI: 12501 NBS:N13:300  Toilet seat, Ergosit, with stabilising buffers, w/cover, raised 50 mm. With raised cross bar hinge	
<b>Colour</b>	White	
<b>Description</b>	Toilet seat, Ergosit, with stabilising buffers, w/cover, raised 50 mm. With raised cross bar hinge	
	  	
<b>Materials</b>	<b>TOILET SEAT:</b> The material is colour ingrained duroplast (UF A 10 = urea formaldehyde) that contains no environmentally hazardous substances. UF plastic comprises approx. 70% urea formaldehyde resin, 30% cellulose, as well as 5% minerals, pigments, lubricants and moisture content. Free formaldehyde is released during the production process. Pressalit toilet seats have been measured after curing (< 0.01 mg/m³) to emit significantly below the permitted EU limit of 0.08 mg/m³.  <b>BUFFERS:</b> EVA (copolymer made of ethylene and vinyl acetate).  <b>MOUNT/HINGES:</b> Stainless steel and thermoplastic.	
<b>Product capacity</b>	This product has been tested statically up to 375 kg. As the European standard for technical aids for disabled, BS/EN 12182:2012, requires a static test with safety factor 1.5, the nominal maximum load for this product is therefore 250 kg.	
<b>Weight</b>	3,5 kg	
<b>Gross weight</b>	3,8 kg	
<b>Weight per crate</b>	0 kg	





# PRESSALIT®

---

**Toilet seat Ergosit with cover, 50 mm raised**  
**With stabilising buffers and special hinges adapted to raised toilet seat.**

---

**Qty per crate** 0

**Qty per pallet** 0

**Cleaning** Use a mild soap solution to clean the toilet seat. Please take care that the seat and hinges are not left damp. For disinfecting the seat, e.g. ethanol 85% can be used. Make sure that the product is wiped down with a soft rag wrung out in clear water and dried off with a soft, dry cloth immediately after disinfection. Avoid contact of the seat and hinges with abrasive, corrosive or chlorine-based cleaners, as these may cause damage or lead to flash rust. Therefore, when cleaning the bowl, make sure that the seat and cover are kept in an upright position until all cleaning agent has been flushed away.

**Warranty** Your Pressalit toilet seat is subject to a 10-year product warranty against presentation of receipt. You should contact the dealer or store where you purchased the seat in the event of a complaint. In the event of complaints about defective materials and/or manufacturing defects and/or functional limitations ascertained within a 10-year period, following proper use and cleaning of the toilet seat (we reserve the right to examine the toilet seat at our factory), Pressalit will deliver spare parts or replace the toilet seat by an identical or similar model. If the color is no longer available, Pressalit reserves its right to replace the toilet seat by a similar model in the color white. Upon replacement or delivery of a spare part, the warranty period will not be extended nor renewed. If replacement is not possible, Pressalit may instead refund the purchase price.

Subject to changes in the product range, technical changes and printing errors.

---

[pressalit.com](http://pressalit.com)

---