




PRESSALIT®

Product Item number	VALUE fold down grab bar, height adjustable R1110000
Color Product description	VVS:615722405 NRF:6191643 RSK:7972412 EAN:5708590176420 HMI: 25413 NBS:N13:436 FDA:ILS/D135866 White  VALUE fold down grab bar, wall-mounted, foldable, height adjustable by 5.9", length 34.1". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately. VALUE fold down grab bar, wall-mounted, foldable, height adjustable by 6", length 34". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, diecast aluminium, plastic. Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.
Adjustability Adjustment height Product capacity	Height adjustable 5,91 inch Maximum load: 265 lbs. When the fold down grab bar is loaded with 265 lbs, each mounting bolt is loaded with 190 lbs.
Materials	Basic material: Steel Surface material: Rilsan End cap: Nylon (PA)

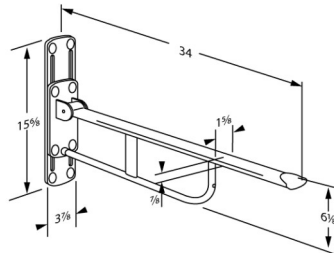

Product	VALUE fold down grab bar, height adjustable
Item number	R1110000
Ergonomy	Provides firm support when you are sitting, standing, or transferring to and from a toilet. Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.
Features	The grab bar can fold up when not in use or to make more space. Height adjustable 6".
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged. For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean) NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.



PRESSALIT®

Product	VALUE fold down grab bar, height adjustable
Item number	R1110000

Mounting in different wall types

	Concrete >B25 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Brick >= Mz12 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Aerated Concrete P4 V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	



Product Item number	VALUE fold down grab bar, height adjustable R1110000
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.

It is the responsibility of the installer to ensure that suitable wall fasteners are used.