




PRESSALIT®

Product Item number	PLUS folding seat, manually height adjustable R5500112
Color	VVS:777531712 NRF:7501063 EAN:5708590328249 HMI: 130162 FDA: INN/D135865 Anthracite gray
Product description	 PLUS folding seat. Manually height adjustable by 9.4" The seat is foldable. For fitting to vertical PLUS wall track (not included). The folding seat is not suitable for use in the shower area. PLUS folding seat. Manually height adjustable by 9.4" Materials: Aluminium, thermoset plastic, plastic. For fitting to vertical track (not included).
Adjustability Adjustment height Product capacity	Height adjustable 9,46 inch Maximum user mass: 584 lbs. When the seat is loaded with 584 lbs, each wall fastener is loaded with 51 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy Features	The folding seat is made out of thermoset plastic, and the mounting and suspension bracket are made from 3/8" anodised aluminium. The end caps are made from ABS plastic. Ergonomically shaped seat, which offers good and comfortable support. The seat is suitable for situations where additional seating is required and space is limited, - e.g. in kitchens, lifts and hallways. The seat can be folded away when not in use, - or to free up space. Stays up by itself when folded away.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

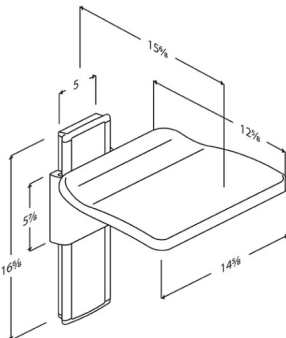

PRESSALIT®

Product	PLUS folding seat, manually height adjustable
Item number	R5500112
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate. DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.



PRESSALIT®

Product	PLUS folding seat, manually height adjustable
Item number	R5500112

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®