




|   |   |
|---|---|
| <p><b>Product</b><br/><b>Item number</b></p>  | <p>PLUS folding seat, manually height adjustable<br/>R5500000</p>   |
| <p><b>Color</b><br/><br/><b>Product description</b></p>                                   | <p>VVS:777531700<br/>NRF:7501062<br/>EAN:5708590328218<br/>HMI: 130162<br/>FDA:INN/D135865<br/>White</p>  <p>PLUS folding seat. Manually height adjustable by 9.4"<br/>The seat is foldable.<br/>For fitting to vertical PLUS wall track (not included).<br/>The folding seat is not suitable for use in the shower area.</p> <p>PLUS folding seat. Manually height adjustable by 9.4"<br/>Materials: Aluminium, thermoset plastic, plastic.<br/>For fitting to vertical track (not included).</p> |
| <p><b>Adjustability</b><br/><b>Adjustment height</b><br/><br/><b>Product capacity</b></p> | <p>Height adjustable<br/>9,46 inch</p> <p>Maximum user mass: 584 lbs.<br/>When the seat is loaded with 584 lbs, each wall fastener is loaded with 51 lbs.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>   |
| <p><b>Materials</b><br/><br/><b>Ergonomy</b><br/><br/><b>Features</b></p>                 | <p>The folding seat is made out of thermoset plastic, and the mounting and suspension bracket are made from 3/8" anodised aluminium. The end caps are made from ABS plastic.</p> <p>Ergonomically shaped seat, which offers good and comfortable support. The seat is suitable for situations where additional seating is required and space is limited, - e.g. in kitchens, lifts and hallways.</p> <p>The seat can be folded away when not in use, - or to free up space. Stays up by itself when folded away.</p>  |
| <p><b>Warranty</b></p>  | <p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>  |

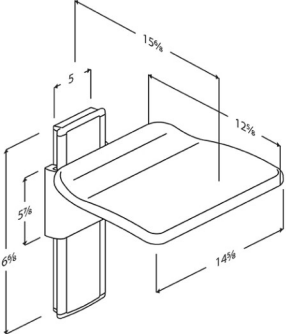

|                    |   |
|--------------------|---|
| <b>Product</b>     | PLUS folding seat, manually height adjustable   |
| <b>Item number</b> | R5500000  |
| <b>Cleaning</b>    | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>Clinimax Difficil-S<br/>Erisan OXY+<br/>Sprit (husholdningsprit)<br/>Suma Bac D10<br/>Virkon<br/>Rodalon Desinfektionsspray<br/>Lysoformin Spezial<br/>CaviCide<br/>Sufanios Premium<br/>Tvättsprit 70%</p> <p><b>NON-DISINFECTANTS:</b><br/>Suma Star free D1<br/>Taski 100<br/>Taski 200<br/>Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b><br/>Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |