




<p>Product Item number</p>	<p>PLUS wash basin bracket, manually sideways adjustable R4560</p>
<p>Colour Product description</p>	<p>VVS:654001300 NRF:6132444 EAN:5708590309408 HMI: 106242 NBS:N13:30,335 FDA:ILW/D142473 Pale grey</p>  <p>PLUS wash basin brackets are available in electrically, pneumatically and manually height adjustable versions. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.</p> <p>PLUS wash basin bracket. Sideways adjustable. Recommendation: Flexible feed and waste system R2064/R2065/R2067. (UK markets R2074/R2075/R2077). Materials: Aluminium and ABS polymer. For fitting to horizontal PLUS wall track (not included).</p>
<p>Adjustability Product capacity</p>	<p>Sideways adjustable</p> <p>Maximum load: 135 kg. Suitable for wash basins weighing no more than 22 kg. When the wash basin is loaded with 135 kg, each wall fastener will be carrying 97 kg. <i>Wash basin brackets need to be able to carry the weight of the wash basin, as well as support the weight of the person who is leaning against the wash basin or who is using the wash basin to rise from seated position.</i> Neither wash basin nor wash basin bracket are intended to carry the full body weight of a person, and the test methods applicable to wash basin brackets reflect this: ISO 17966 prescribes a static test on the front edge of the wash basin at 75% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 135 kg and the wash basin weighs 22 kg, it is – according to ISO 17966 – the equivalent of a user weighing 151 kg leaning on the front edge of the wash basin.</p>
<p>Materials Ergonomy</p>	<p>The wash basin bracket is primarily made of 10-micrometre anodized aluminium. The ends of the bracket cover are capped with ABS polymer.</p> <p>The wash basin bracket suits all types of bolt-mounted wash basin. Its rounded, sleek design greatly facilitates cleaning.</p>



Product	PLUS wash basin bracket, manually sideways adjustable
Item number	R4560
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>PLUS & SELECT wash basin brackets The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
Declaration of conformity	<p>CE UK MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices. For the UK: - Medical Device Regulation (MDR) 2002 on medical devices.</p>

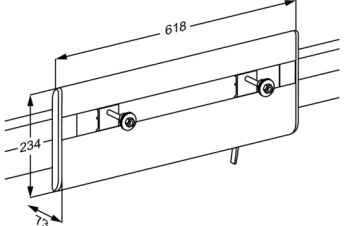

PRESSALIT®

Product	PLUS wash basin bracket, manually sideways adjustable
Item number	R4560

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®