




<p><b>Product</b></p> <p><b>Item number</b></p>	<p>PLUS fold down grab bar with integrated counter-balance, 33.5", right hand operated</p> <p>R381485277</p>
<p><b>Color</b></p> <p><b>Product description</b></p>	<p>NRF:6191154  EAN:5708590368801  HMI: 123635  FDA: ILS/D135866</p> <p>Lime green</p>  <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions. Suitable for use with elongated toilet pans or floor standing toilets. <i>Replaces discontinued model R3285.</i></p> <p><i>PLUS fold down grab bar. Projection 34.5", ergonomically shaped to provide a clear sense of where the fold down grab bar ends. adjustable horizontally. Detachable and foldable.</i></p> <p><i>Materials: The arms are constructed of 3/8" anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core For mounting on horizontal wall track (not included).</i></p>
<p><b>Adjustability</b></p> <p><b>Adjustment height</b></p> <p><b>Product capacity</b></p>	<p>Horizontally adjustable</p> <p>9,46 inch</p> <p>Maximum load: 330 lbs.  When the fold down grab bar is loaded with 330 lbs, each wall fastener in the horizontal wall track is loaded with 16 kg.</p> <p><i>PLUS fold down grab bars have been tested statically with 1540 lbs (close to the wall), 935 lbs (in the middle) and 495 lbs (the outermost point).</i></p>
<p><b>Materials</b></p>	<p>The grab bars are primarily constructed in 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core.</p> <p>The hinge arm is also made of soft SEBS.</p> <p>Internal components and fasteners are primarily made of stainless steel.</p> <p>The hinge barrel is capped with ABS polymer.</p>

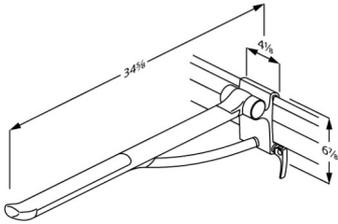

# PRESSALIT®

<b>Product</b>	PLUS fold down grab bar with integrated counter-balance, 33.5", right hand operated
<b>Item number</b>	R381485277
<b>Ergonomy</b>	The oval shape gives an easy and comfortable surface to lean and push against. The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends. Soft and slip-resistant material that has a warm feel.
<b>Features</b>	Folds against the wall. The underside is shaped for the mounting of various accessories.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



<b>Product</b>	PLUS fold down grab bar with integrated counter-balance, 33.5", right hand operated
<b>Item number</b>	R381485277
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



<b>Product</b>	PLUS fold down grab bar with integrated counter-balance, 33.5", right hand operated
<b>Item number</b>	R381485277

**Intended use**

The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.

**Operating environment**

For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals.

For institutional use, the expected life time is 10 years.

*Life cycle test performed according to ISO 17966:2016.*

For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".