




<b>Product</b> <b>Item number</b>	VALUE IV shower seat, fixed height R1620000112
<b>Color</b>	VVS:777525312 NRF:6314097 RSK:8162781 EAN:5708590377711 HMI: 124483 DB: 2296038 FDA:INN/D135865 Anthracite gray
<b>Product description</b>	 <p>VALUE IV folding shower seat and backrest create a simple solution for users who need to sit down while showering.</p> <p>VALUE IV shower seat. Fixed height. The seat is foldable.          Materials: Diecast aluminium, Rilsan-coated stainless steel and vacuum-moulded polyethylene (PE).          Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p>
<b>Adjustability</b> <b>Product capacity</b>	Fixed mounting Maximum user mass: 298 lbs. When the seat is loaded with 298 lbs, each wall fastener is loaded with 276 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
<b>Materials</b>	The wall profile is made of diecast aluminium, and the frame is made of Rilsan-coated stainless steel. Rilsan provides good protection against corrosion and is comfortable to touch and easy to clean. The seat is made of vacuum-moulded polyethylene (PE).
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

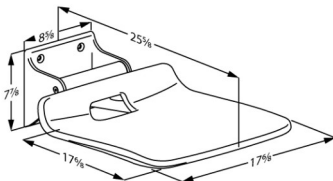

# PRESSALIT®

<b>Product</b>	VALUE IV shower seat, fixed height
<b>Item number</b>	R1620000112
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b></p> <ul style="list-style-type: none"><li>Clinimax Difficil-S</li><li>Erisan OXY+</li><li>Sprit (husholdningsprit)</li><li>Suma Bac D10</li><li>Virkon</li><li>Rodalon Desinfektionsspray</li><li>Lysoformin Spezial</li><li>CaviCide</li><li>Sufanios Premium</li><li>Tvättspnit 70%</li><li>WipeClean UniversalCleaning</li><li>Chlorine Disinfection Wipes</li><li>Klorilli (Kiilto Clean)</li></ul> <p><b>NON-DISINFECTANTS:</b></p> <ul style="list-style-type: none"><li>Suma Star free D1</li><li>Taski 100</li><li>Taski 200</li><li>Milizid Citro</li></ul> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	<p>This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.</p>



# PRESSALIT®

<b>Product</b>	VALUE IV shower seat, fixed height
<b>Item number</b>	R1620000112

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®