



# PRESSALIT®

<b>Product</b> <b>Item number</b>	VALUE fold down grab bar floor mounted on a column R1161000
<b>Color</b>	VVS:615722425 NRF:6190517 RSK:8980351 EAN:5708590176697 HMI: 25416 NBS:N13:436 FDA:ILS/D135866 White
<b>Product description</b>	<p>VALUE fold down grab bar, floor-mounted on free-standing column, length 37.6". Height adjustable by 26.4". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. fixing bolts for solid floorl. Plugs depend on floor type and must be ordered separately.</p> <p>VALUE fold down grab bar, floor-mounted on free-standing column, length 37.6". Height adjustable by 26.4". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, pressure-cast aluminium, plastic. Incl. fixing bolts for solid floor. Plugs depend on floorl type and must be ordered separately.</p>
<b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b>	Height adjustable 26,4 inch Max. vertical load: 265 lbs.
<b>Materials</b> <b>Ergonomy</b>	Core material: Steel Surface material: Rilsan End cap: Nylon (PA) Provides firm support when you are sitting, standing, or transferring to and from a toilet. A wide support surface to hold on to. The rounded shape provides a firm and comfortable grip.



# PRESSALIT®

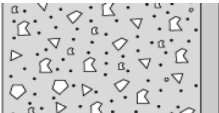
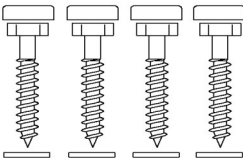

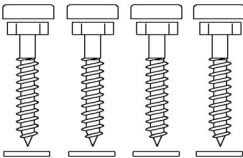
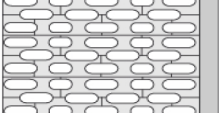
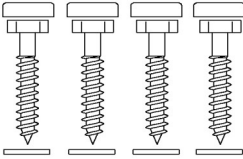
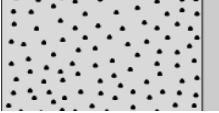
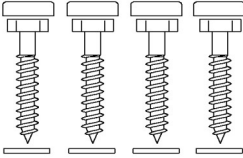
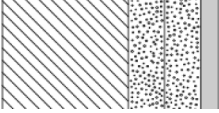
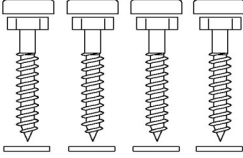
<b>Product</b>	VALUE fold down grab bar floor mounted on a column
<b>Item number</b>	R1161000
<b>Features</b>	The space saving grab bars can fold up when not in use. Height adjustable.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean)</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>



PRESSALIT®

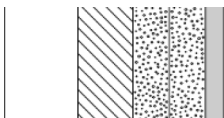
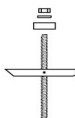

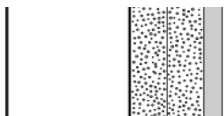
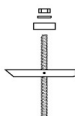

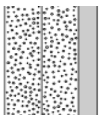
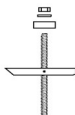

<b>Product</b>	VALUE fold down grab bar floor mounted on a column
<b>Item number</b>	R1161000

## Mounting in different wall types

	<p><b>Concrete &gt;B25</b></p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p><b>Solid Brick &gt;= Mz12</b></p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p><b>Solid Sand Lime Brick &gt;= Hlz12</b></p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p><b>Aerated Concrete P4</b></p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b></p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	



# PRESSALIT®

Product Item number	VALUE fold down grab bar floor mounted on a column R1161000		
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b> Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b> Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	<b>Drywall (2 x 12,5 mm) water-resistant</b> Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.

[pressalit.com](http://pressalit.com)

Pressalit A/S  
 Pressalitvej 1  
 DK-8680 Ry  
 Tel.: +45 8788 8788  
 Fax: +45 8788 8789  
 E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®