




PRESSALIT®

Product Item number	TMA 4 arm support for mounting on toilet (no call button) R1196112
Colour	VVS:615726112 NRF:1617307 EAN:5708590397160 Anthracite grey 
Product description	TMA arm supports comes in various widths, and models with or without toilet seat. TMA 4 arm supports for mounting on toilet. Anthracite grey handles. Incl. two paper holders. Excl. hinges. Bottom mounting. Compatible with Gustavsberg Nautic 1546 toilet Recommended with toiletseats: Pressalit Sign 553011-DN1999 Pressalit Sign 553112-DN1999 Pressalit Sign 554011-DN1999 Pressalit Sign 554112-DN1999
Adjustability Product capacity	Fixed mounting Maximum user weight (MUM): 200 kg. The armsupports are tested according to ISO/FDIS 17966:2025
Materials Ergonomy Features	Stainless steel Aluminium Polyurethane The arm supports provide firm support and the handles are comfortable for the hands and arms. The TMA 4 features individually foldable arm supports, allowing flexible adjustment and easy access tailored to different user needs and care scenarios. The arm supports provide firm and reliable stability, while the ergonomically designed handles ensure a high level of comfort for both hands and arms during use.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

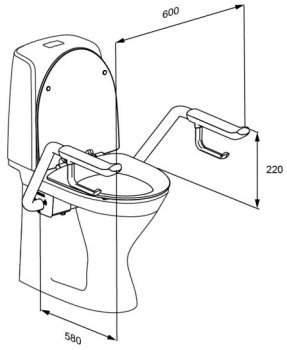

PRESSALIT®

Product	TMA 4 arm support for mounting on toilet (no call button)
Item number	R1196112
Cleaning	The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.

pressalit.com
