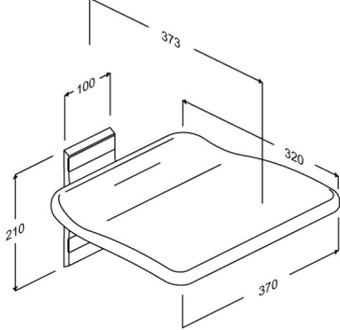




<p>Product Item number</p>	<p>PLUS folding seat, fixed height R5510112</p>
<p>Colour Product description</p>	<p>VVS:777531012 NRF:7501002 RSK:8565072 EAN:5708590321332 HMI: 130159 NBS:N13:448 FDA:INN/D135865 Anthracite grey</p>  <p>The seat is suitable for situations where additional seating is required and space is limited, - e.g. in kitchens, lifts and corridors or on the landings. The folding seat is not suitable for use in the shower area.</p> <p>PLUS folding seat. For fixed mounting. Materials: Aluminium, thermoset plastic, plastic. Excl. mounting kit.</p>
<p>Adjustability Product capacity</p>	<p>Fixed mounting</p> <p>Maximum user mass: 435 kg. When the seat is loaded with 435 kg, each wall fastener is loaded with 199 kg.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>
<p>Materials Ergonomy Features</p>	<p>The folding seat is made out of thermoset plastic, and the mounting and suspension bracket are made from 10 my anodised aluminium. The end caps are made from ABS plastic.</p> <p>Ergonomically shaped seat, which offers good and comfortable support.</p> <p>The seat can be folded away when not in use, - or to free up space. Stays up of itself when folded away.</p>
<p>Warranty</p>	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>

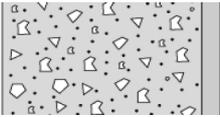
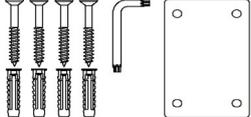
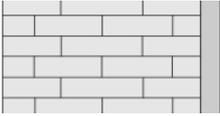
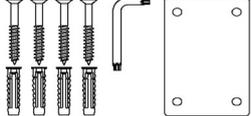
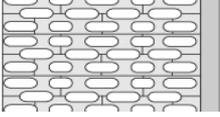
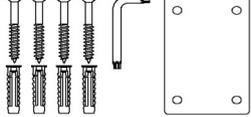
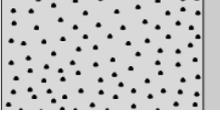
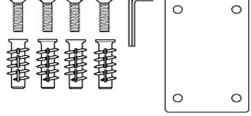
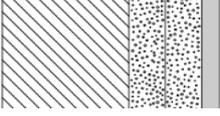
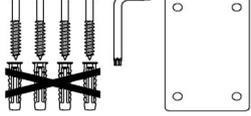


Product	PLUS folding seat, fixed height
Item number	R5510112
Cleaning	<p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

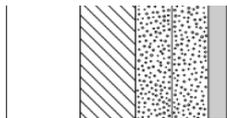
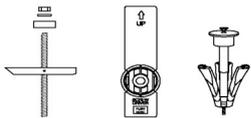
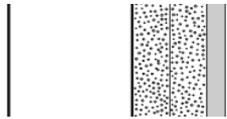
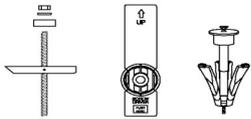
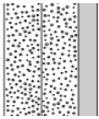


Product	PLUS folding seat, fixed height
Item number	R5510112

Installation in different wall types

	Concrete >B25 V8604 - Mounting kit, 4 countersunk screws Ø8 x 70, plugs SX10. To be ordered separately.	
	Solid Brick >= Mz12 V8604 - Mounting kit, 4 countersunk screws Ø8 x 70, plugs SX10. To be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8604 - Mounting kit, 4 countersunk screws Ø8 x 70, plugs SX10. To be ordered separately.	
	Aerated Concrete P4 V8654 - Mounting kit, 4 countersunk bolts M8 x 35, plugs FTP M 8. To be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8604 - Mounting kit, 4 countersunk screws Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.	



Product Item number	PLUS folding seat, fixed height R5510112
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. 
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. 
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. 