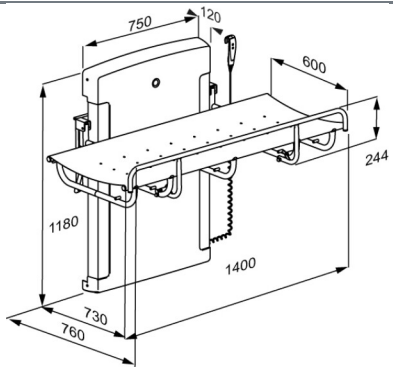




PRESSALIT®

Product Item number	SCT 1000 shower change table, electrically height adjustable R8583114000
Colour	VVS: 777543640 EAN: 5708590377636 HMI: 140128 White
Product description	 <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>SCT 1000 showerchange table, powered, height adjustable 400 mm; from 480 mm - 880 mm, with wired hand control Foldable. Incl. safety rail. Outer measurements 760 mm x 1400 mm. Excl. mounting kit.</p>
Adjustability Adjustment height Product capacity Speed upwards Noise level	Height adjustable 400 mm Maximum user weight (MUM): 200 kg. Safe working load (SWL): 200 kg. When the table is loaded with 200 kg, each mounting bolt supports 27 kg. <i>This product has been tested statically with a load of 500 kg as the IEC 60601-1 standard requires testing with safety factor 2.5.</i> Speed up mm: 11 mm/s Speed down mm: 9 mm/s 50 dB
Materials	The table is made from practical and durable PVC canvas which is phthalate free. The frame is powder-coated stainless steel.



Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8583114000
Ergonomy	Carers can easily set the height of the table, enabling them to obtain the correct working position. This ensures safe showering and changing for the user as well as minimizing the risk of back injury for the carer.
Features	Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall. The integrated gas cylinders enable the table to be folded out of the way when not in use to maximise floor space in front of the unit. Splash curtain (optional) for the shower and changing table protects the carer against splashing.
Mains connection (V)	230
Actuator (motor)	LA40 Careline (Linak), 4000 N, IPX6
Control box	CA63 Careline (Linak), IPX6
Wired remote control	HB30 (Linak), IP66
Power consumption (maximum load)	Power consumption by maximum load, watt: 420 W Power consumption by standby: 0,5 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



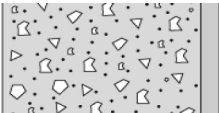
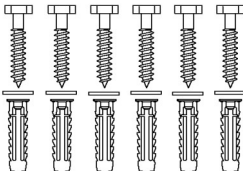

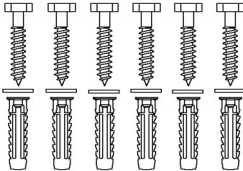
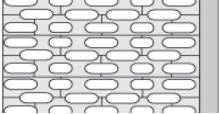
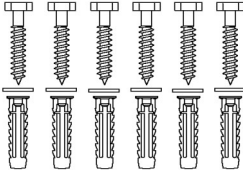
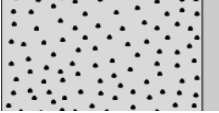
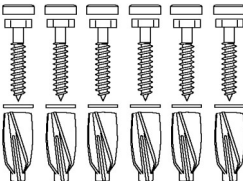
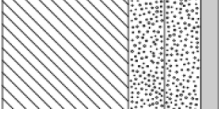
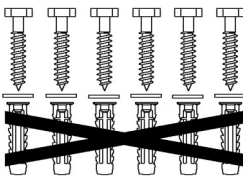
Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8583114000
Cleaning	<p>SCT 1000 & 2000</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes OxiGenesis from ENVIROX</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>



PRESSALIT®

Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8583114000

Installation in different wall types

	<p>Concrete >B25</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	



PRESSALIT®

Product Item number	SCT 1000 shower change table, electrically height adjustable R8583114000
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant</p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®