




PRESSALIT®

Product Item number	PLUS replacement shower seat 450, manually height adjustable R7134000
Color	VVS:777514800 NRF:6317076 EAN:5708590359151 FDA:INN/D135865 White
Product description	 PLUS shower seats are offered in a number of colours. They are available as both fixed height, electrically adjustable and manually adjustable, with and without aperture, and with an option of back and armrests. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks. PLUS shower seat 450, with back- and armrests. Manually height adjustable by 7.7", without the use of tools. Armrests are independently height adjustable. The seat and the armrests are foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS incorporating a glass-reinforced nylon core. For mounting on vertical wall track (not included).
Adjustability Adjustment height Product capacity	Height adjustable 7,68 inch Maximum user mass (seat): 440 lbs. Maximum load for armrests is 55 kg. When the seat is loaded with 440 lbs, each wall fastener is loaded with 40 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials	The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.

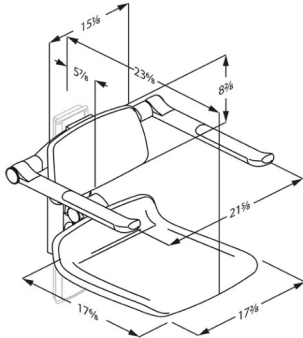

Product Item number	PLUS replacement shower seat 450, manually height adjustable R7134000
Ergonomy Features	Ergonomically designed seat - contact surface are 12.2" x 14.4". Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against. The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximise floor space.
Warranty Cleaning	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Erisan OXY+ Tvättspit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH) NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa) Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water. The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate. DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.



PRESSALIT®

Product	PLUS replacement shower seat 450, manually height adjustable
Item number	R7134000

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®

8. juli 2026