




|   |   |
|---|---|
| <p><b>Product</b><br/><b>Item number</b></p>                    | <p>PLUS backrest, fixed height<br/>R6304277</p>   |
| <p><b>Color</b><br/><br/><b>Product description</b></p>         | <p>NRF:6317029<br/>EAN:5708590365343<br/>HMI: 124221<br/>FDA:INN/D135865</p> <p>Lime green</p>  <p>PLUS backrests are offered in a number of colours and versions, they are available as both fixed height and adjustable. They are suitable for wall-mounting behind a toilet without cistern or mounting on one of Pressalit's toilet lifters.</p> <p>PLUS backrest, 11.5" x 7.5", soft feel.<br/>Fixed height, for mounting directly on the wall behind a wall-mounted toilet without cistern.<br/>Materials: Aluminium and polyurethane foam with a plywood core.<br/>Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> |
| <p><b>Adjustability</b><br/><b>Product capacity</b></p>         | <p>Fixed mounting</p> <p>Maximum load: 110 lbs.<br/>Maximum user weight: 220 lbs.<br/>When the backrest is loaded with 110 lbs, each wall fastener will be carrying 130 lbs.<br/><i>Backrests need to be able to support the weight of the person who is leaning against it.</i><br/><b><i>The backrest is not intended to carry the full body weight of a person, and the test methods applicable to backrests reflect this: ISO 17966 prescribes a static test at 50% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 110 lbs, it is – according to ISO 17966 – the equivalent of a user weighing 220 lbs leaning on the backrest.</i></b></p>   |
| <p><b>Materials</b><br/><b>Ergonomy</b><br/><b>Features</b></p> | <p>The backrest is made of polyurethane (PU) foam which is covered with a moulded coated surface that is protected using a clear PU lacquer.</p> <p>Soft and slip-resistant material that has a warm feel.</p> <p>For toilet without cistern</p>  |
| <p><b>Warranty</b></p>  | <p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>  |

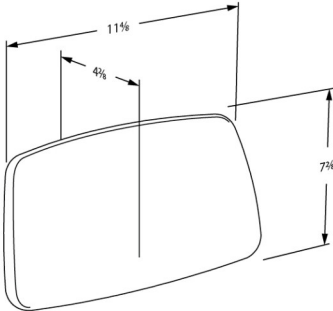

|                    |   |
|--------------------|---|
| <b>Product</b>     | PLUS backrest, fixed height   |
| <b>Item number</b> | R6304277  |
| <b>Cleaning</b>    | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>Erisan OXY+<br/>Tvättsprit 70%<br/>CaviCide surface disinfectant decontaminant cleaner<br/>Rodalon Desinfektionsspray<br/>Terralin protect (SCHÜLKE+)<br/>Disonzon Plus (BODE Chemie GmbH)</p> <p><b>NON-DISINFECTANTS:</b><br/>Wetrok Caledor (Wetrok)<br/>Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b><br/>Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |