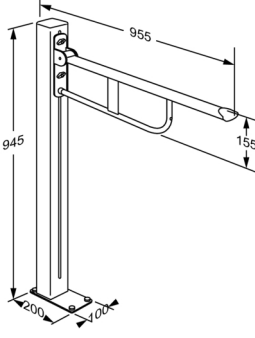




Product	VALUE support arm on column
Item number	R1161000
Colour	VVS:615722425 NRF:6190517 RSK:8980351 EAN:5708590176697 HMI: 25416 NBS:N13:436 FDA:ILS/D135866 White
Product description	 VALUE support arm, floor-mounted on free-standing column, projection 955 mm. Height adjustable 670 mm. The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. fixing bolts for solid floor. Plugs depend on floor type and must be ordered separately. VALUE support arm, floor-mounted on free-standing column, projection 955 mm. Height adjustable 670 mm. The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, pressure-cast aluminium, plastic. Incl. fixing bolts for solid floor. Plugs depend on floor type and must be ordered separately.
Adjustability	Height adjustable
Adjustment height	670 mm
Product capacity	Max. vertical load: 120 kg.
Materials	Basic material: Steel Surface material: Rilsan End cap: Nylon (PA)



PRESSALIT®

Product	VALUE support arm on column
Item number	R1161000
Ergonomy	Provides firm support when you are sitting, standing, or transferring to and from a toilet. Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.
Features	The support arm can fold up when not in use or to make more space. Height adjustable.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

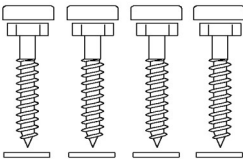
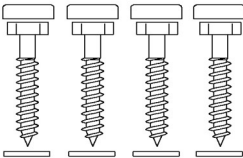
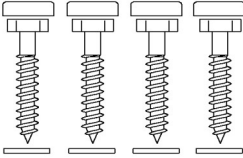
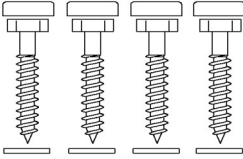
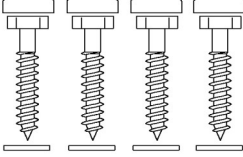


Product	VALUE support arm on column
Item number	R1161000
Cleaning	VALUE support arms The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged. For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean) NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.



Product	VALUE support arm on column
Item number	R1161000

Installation in different wall types

	Concrete >B25 V8596 - 4 coach screws Ø10 x 70, with white screw covers are included. Plugs to be ordered separately.	
	Solid Brick >= Mz12 V8596 - 4 coach screws Ø10 x 70, with white screw covers are included. Plugs to be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8596 - 4 coach screws Ø10 x 70, with white screw covers are included. Plugs to be ordered separately.	
	Aerated Concrete P4 V8596 - 4 coach screws Ø10 x 70, with white screw covers are included. Plugs to be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8596 - 4 coach screws Ø10 x 70, with white screw covers are included. Plugs to be ordered separately.	



PRESSALIT®

Product Item number	VALUE support arm on column R1161000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®