




<p><b>Product</b></p> <p><b>Item number</b></p>	<p>PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated</p> <p>R372470000</p>
<p><b>Color</b></p> <p><b>Product description</b></p>	<p>VVS:615877400            NRF:6191131            EAN:5708590368573            HMI: 138181            FDA:ILS/D135866            White</p>  <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions.</p> <p>PLUS fold down grab bar, counter-balanced, 27.6".            Height adjustable 5.9". Foldable. Right handed.            Only for use with TL3/TL4 toilet lifter.            Mounting profile included.            Ergonomically shaped to provide a clear sense of where the support arm ends.            Materials: The arms are constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core.</p>
<p><b>Adjustability</b></p> <p><b>Adjustment height</b></p> <p><b>Product capacity</b></p>	<p>Height adjustable</p> <p>5,91 inch</p> <p>Maximum load: 330 lbs.</p>
<p><b>Materials</b></p> <p><b>Ergonomy</b></p>	<p>The grab bars are primarily constructed from 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core.            The hinge arm is also made of soft SEBS.            Internal components and fasteners are primarily made of stainless steel.            The hinge barrel is capped with ABS polymer.</p> <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p>



<b>Product</b>	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated
<b>Item number</b>	R372470000
<b>Features</b>	The underside is shaped for the mounting of various accessories.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>

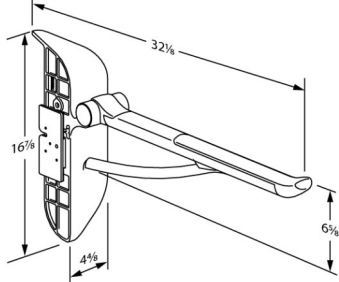

<b>Product</b>	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated
<b>Item number</b>	R372470000

**Intended use**

The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.

**Operating environment**

For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals.

For institutional use, the expected life time is 10 years.

*Life cycle test performed according to ISO 17966:2016.*

For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".