




PRESSALIT®

Product Item number	MSCT 4 shower change trolley, manually height adjustable R9504721000
Color	VVS: 777542630 EAN: 5708590389851 HMI: 147809 FDA: INW/D546021 White
Product description	 Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user. MSCT 4 shower change trolley, 27.6" x 82.7", four fold-down sides. Manually height adjustable with foot pedal (hydraulic). 19.7" lifting range, from 21.7" - 41.3" Mattress to be ordered separately.
Adjustability Adjustment height Product capacity	Height adjustable 19,7 inch Maximum user weight (MUM): 441 lbs. Maximum load (SWL): 551 lbs. <i>The maximum load is to be understood as the total safe working load of the product, including both the weight of the user, and the weight of the mattress</i> <i>This product has been tested statically with a load of 826.5 lbs, as the ISO 17966 standard requires testing with a safety factor 1.5.</i>
Materials	Plastics, Compact Density Fiberboards (CDF) and powder-lacquered stainless steel.

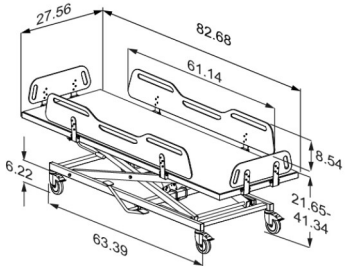

Product Item number	MSCT 4 shower change trolley, manually height adjustable R9504721000
Ergonomy	• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user. • One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table. • The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training. • Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.
Features	• The trolley can be moved from room to room for easy transport and multiple users. • Intuitive, easy-to-use drop down rails with one hand operation.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



PRESSALIT®

Product	MSCT 4 shower change trolley, manually height adjustable
Item number	R9504721000
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.

pressalit.com

PRESSALIT®