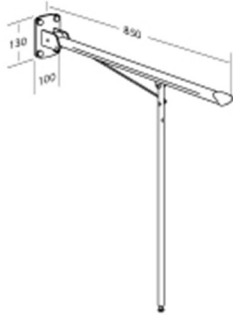




<b>Product</b> <b>Item number</b>	VALUE support arm with leg, fixed height R1101000
<b>Colour</b> <b>Product description</b>	VVS:615722400 NRF:6191641 EAN:5708590176413 HMI: 25412 NBS:N13:436 FDA: ILS/D135866 White  VALUE support arm with leg, fixed, foldable, projection 850 mm, height 850 mm. The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately. 
<b>Adjustability</b> <b>Product capacity</b>	Fixed mounting Maximum load: 135 kg. When the support arm is loaded with 135 kg, each wall fastener is loaded with 81 kg.
<b>Materials</b> <b>Ergonomy</b> <b>Features</b>	The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA). Provides firm support when you are sitting, standing, or transferring to and from a toilet.  Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.  The support arm can fold up when not in use or to make more space.  With leg for extra support.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.





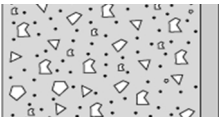
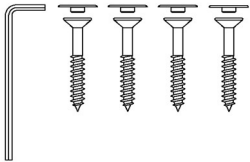

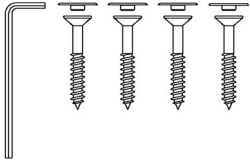
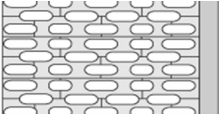
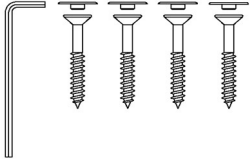
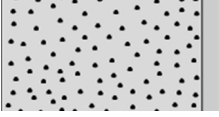
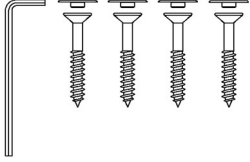
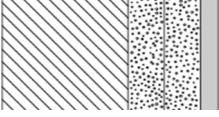
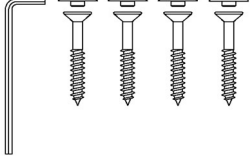
<b>Product</b>	VALUE support arm with leg, fixed height
<b>Item number</b>	R1101000
<b>Cleaning</b>	<p>VALUE support arms. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged. For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Klorilli</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	<p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007.</p>





<b>Product</b>	VALUE support arm with leg, fixed height
<b>Item number</b>	R1101000

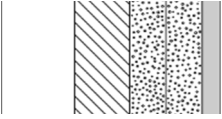
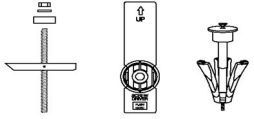
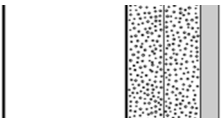
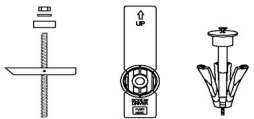
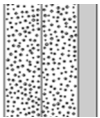

## Installation in different wall types

	<b>Concrete &gt;B25</b> V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	<b>Solid Brick &gt;= Mz12</b> V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	<b>Solid Sand Lime Brick &gt;= Hlz12</b> V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	<b>Aerated Concrete P4</b> V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b> V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.	





# PRESSALIT®

<b>Product</b> <b>Item number</b>	VALUE support arm with leg, fixed height R1101000
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b> Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. 
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</b> Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. 
	<b>Drywall (2 x 12,5 mm) water-resistant</b> Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. 

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.



[www.pressalit.com](http://www.pressalit.com)

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®