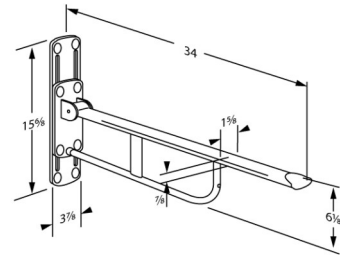




PRESSALIT®

Product	VALUE fold down grab bar, height adjustable
Item number	R1110000
Color	<p>VVS:615722405 NRF:6191643 RSK:7972412 EAN:5708590176420 HMI: 25413 NBS:N13:436 FDA:ILS/D135866</p> <p>White</p>
Product description	<p>VALUE fold down grab bar, wall-mounted, foldable, height adjustable by 5.9", length 34.1". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> <p>VALUE fold down grab bar, wall-mounted, foldable, height adjustable by 6", length 34". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, diecast aluminium, plastic. Incl. mounting bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p>
Adjustability	Height adjustable
Adjustment height	5,91 inch
Product capacity	<p>Maximum load: 265 lbs. When the fold down grab bar is loaded with 265 lbs, each mounting bolt is loaded with 190 lbs.</p>
Materials	<p>Basic material: Steel</p> <p>Surface material: Rilsan</p> <p>End cap: Nylon (PA)</p>





PRESSALIT®

Product	VALUE fold down grab bar, height adjustable
Item number	R1110000
Ergonomy	<p>Provides firm support when you are sitting, standing, or transferring to and from a toilet.</p> <p>Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.</p>
Features	<p>The grab bar can fold up when not in use or to make more space.</p> <p>Height adjustable 6".</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



PRESSALIT®

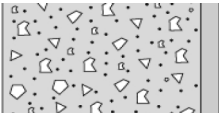
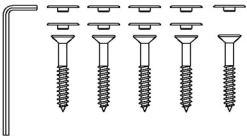

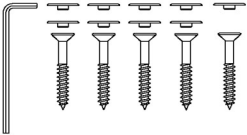
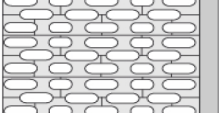
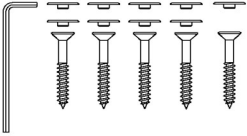
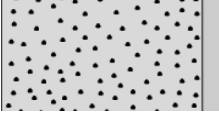
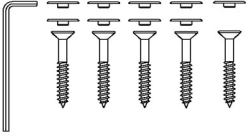
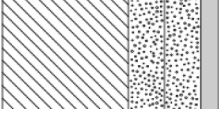
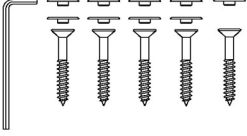
Product	VALUE fold down grab bar, height adjustable
Item number	R1110000
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean)</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
Declaration of conformity	<p>This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.</p>



PRESSALIT®

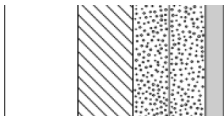
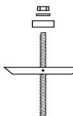

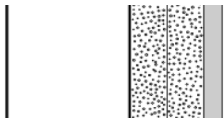
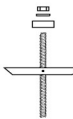

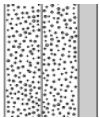

Product	VALUE fold down grab bar, height adjustable
Item number	R1110000

Mounting in different wall types

	<p>Concrete >B25</p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8592 - 5 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	



PRESSALIT®

Product Item number	VALUE fold down grab bar, height adjustable R1110000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com