



PRESSALIT®

Product Item number	VALUE fold down grab bar floor mounted on a column R1161000
Color	VVS:615722425 NRF:6190517 RSK:8980351 EAN:5708590176697 HMI: 25416 NBS:N13:436 White
Product description	<p>VALUE fold down grab bar, floor-mounted on free-standing column, length 37.6". Height adjustable by 26.4". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. fixing bolts for solid floorl. Plugs depend on floor type and must be ordered separately.</p> <p>VALUE fold down grab bar, floor-mounted on free-standing column, length 37.6". Height adjustable by 26.4". The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, pressure-cast aluminium, plastic. Incl. fixing bolts for solid floor. Plugs depend on floorl type and must be ordered separately.</p>
Adjustability Adjustment height Product capacity	Height adjustable 26,4 inch Max. vertical load: 265 lbs.
Materials Ergonomy Features	Core material: Steel Surface material: Rilsan End cap: Nylon (PA) Provides firm support when you are sitting, standing, or transferring to and from a toilet. A wide support surface to hold on to. The rounded shape provides a firm and comfortable grip. The space saving grab bars can fold up when not in use. Height adjustable.



PRESSALIT®

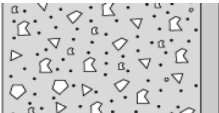
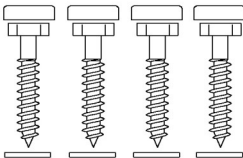

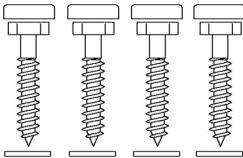
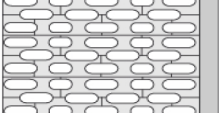
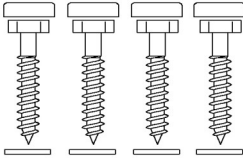
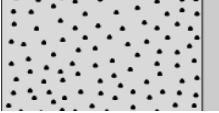
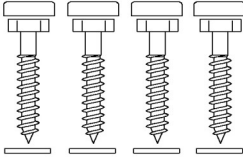
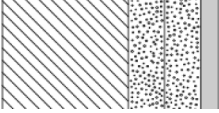
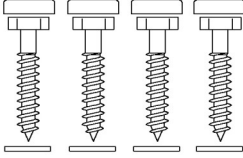
Product	VALUE fold down grab bar floor mounted on a column
Item number	R1161000
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean)</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
Declaration of conformity	This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.



PRESSALIT®

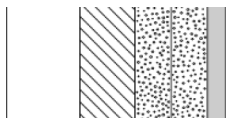
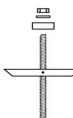


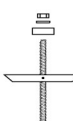

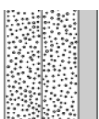
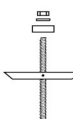

Product	VALUE fold down grab bar floor mounted on a column
Item number	R1161000

Mounting in different wall types

	<p>Concrete >B25</p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8596 - 4 lag bolts Ø10 x 70, with white covers are included. Plugs to be ordered separately.</p>	



PRESSALIT®

Product Item number	VALUE fold down grab bar floor mounted on a column R1161000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®