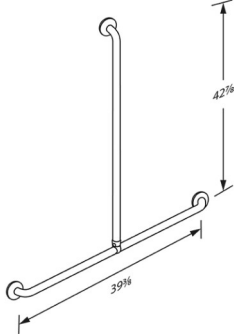




|   |  |
|---|--|
| <p><b>Product</b><br/><b>Item number</b></p>            | <p>PLUS grab bar combination for shower, 39.4" x 42.9"<br/>RT144000</p>  |
| <p><b>Color</b><br/><br/><b>Product description</b></p> | <p>VVS:777569700<br/>NRF:6311235<br/>RSK:8180962<br/>EAN:5708590247649<br/>HMI: 42068<br/>NBS:N13:436</p> <p>White</p>  <p>PLUS grab bars are high-quality products and they can be combined in countless ways and they provide the final touch to furnishing a functional and flexible bathroom. A reliable solution with many different options, which can easily lay the grounds for increased independence.</p> <p>PLUS grab bar combination for shower, A: 39.4" x B: 42.9".<br/>Compatible with RT600 shower holder.<br/>Material: Powder-coated aluminium and plastic.<br/>Excl. mounting kit.</p> |
| <p><b>Product capacity</b></p>                          | <p>Maximum load: 440 lbs.<br/>When the grab bar is loaded with 440 lbs, each wall fastener is loaded with 123 lbs.<br/><i>Grab bars are tested with a 1.5 safety factor in accordance with the requirements in ISO 17966:2016.</i></p>   |
| <p><b>Materials</b></p>                                 | <p>Core material: Aluminium<br/>Surface material: Polyester lacquer<br/>Cover cap: Cycloy plastic<br/>Strap: Nylon (PA)</p>  |
| <p><b>Warranty</b></p>                                  | <p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>   |

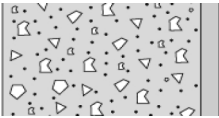
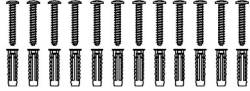

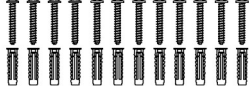
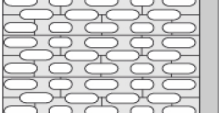
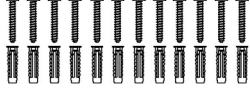
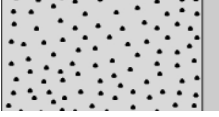

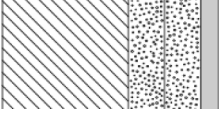
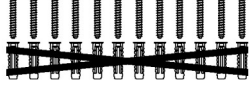


|                                  |   |
|----------------------------------|---|
| <b>Product</b>                   | PLUS grab bar combination for shower, 39.4" x 42.9"   |
| <b>Item number</b>               | RT144000  |
| <b>Cleaning</b>                  | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b><br/>Clinimax Difficil-S<br/>Erisan OXY+<br/>Sprit (husholdningsprit)<br/>Suma Bac D10<br/>Virkon<br/>Rodalon Desinfektionsspray<br/>Lysoformin Spezial<br/>CaviCide<br/>Sufanios Premium<br/>Tvättsprit 70%</p> <p><b>NON-DISINFECTANTS:</b><br/>Suma Star free D1<br/>Taski 100<br/>Taski 200<br/>Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b><br/>Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |
| <b>Declaration of conformity</b> | <p>This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.</p>   |

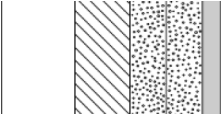
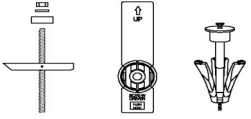
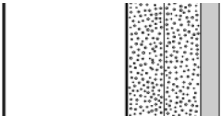
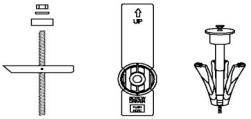
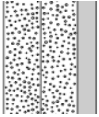



|                    |   |
|--------------------|---|
| <b>Product</b>     | PLUS grab bar combination for shower, 39.4" x 42.9" |
| <b>Item number</b> | RT144000  |

## Mounting in different wall types

|   |  |   |
|---|--|---|
|    | <b>Concrete &gt;B25</b><br>V8617 - Mounting kit, 12 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.   |    |
|   | <b>Solid Brick &gt;= Mz12</b><br>V8617 - Mounting kit, 12 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.   |   |
|  | <b>Solid Sand Lime Brick &gt;= Hlz12</b><br>V8617 - Mounting kit, 12 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.  |  |
|  | <b>Aerated Concrete P4</b><br>V8667 - Mounting kit, 12 bolts M5 x 50 - plugs FTP K 6, for aerated concrete. To be ordered separately.  |  |
|  | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b><br>V8617 - Mounting kit, 12 bolts Ø5 x 50 - plugs SX8. NB!<br>Plugs are not needed.<br>To be ordered separately. |  |



|   |  |   |
|---|--|---|
| <b>Product</b><br><b>Item number</b>  | PLUS grab bar combination for shower, 39.4" x 42.9"<br>RT144000  |   |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b><br>Wall fasteners (12 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or Wingt anchors or similar.           |    |
|    | <b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b><br>Wall fasteners (12 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or Wingt anchors or similar. |    |
|  | <b>Drywall (2 x 12,5 mm) water-resistant</b><br>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.   |  |