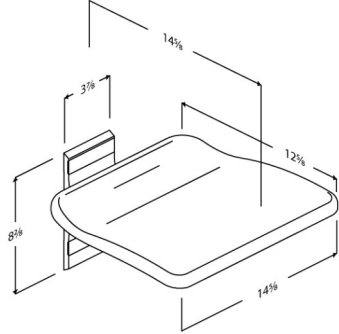




PRESSALIT®

Product Item number	PLUS folding seat, fixed height R5510000
Color	VVS:777531000 NRF:7501001 RSK:7973939 EAN:5708590306339 HMI: 130159 NBS:N13:448 FDA:INN/D135865 White
Product description	 <p>The seat is suitable for situations where additional seating is required and space is limited, - e.g. in kitchens, lifts and hallways. The folding seat is not suitable for use in the shower area.</p> <p>PLUS folding seat. For fixed mounting. Materials: Aluminium, thermoset plastic, plastic. Excl. mounting kit.</p>
Adjustability Product capacity	Fixed mounting Maximum user mass: 960 lbs. When the seat is loaded with 960 lbs, each wall fastener is loaded with 436 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy Features	The folding seat is made out of thermoset plastic, and the mounting and suspension bracket are made from 3/8" anodised aluminium. The end caps are made from ABS plastic. Ergonomically shaped seat, which offers good and comfortable support. The seat can be folded away when not in use, - or to free up space. Stays up by itself when folded away.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



PRESSALIT®

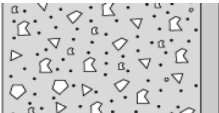
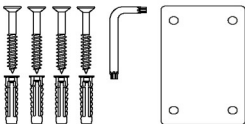

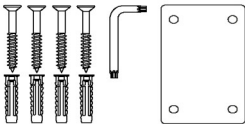
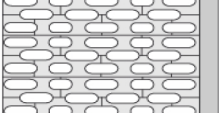
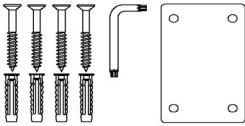
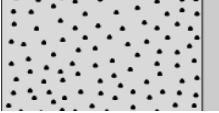
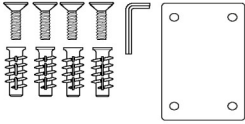
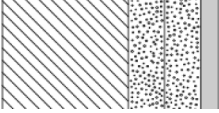
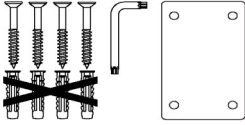
Product	PLUS folding seat, fixed height
Item number	R5510000
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>
Declaration of conformity	<p>This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.</p>



PRESSALIT®

Product	PLUS folding seat, fixed height
Item number	R5510000

Mounting in different wall types

	<p>Concrete >B25</p> <p>V8604 - Mounting kit, 4 countersunk bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Brick >= Mz12</p> <p>V8604 - Mounting kit, 4 countersunk bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8604 - Mounting kit, 4 countersunk bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Aerated Concrete P4</p> <p>V8654 - Mounting kit, 4 countersunk bolts M8 x 35, plugs FTP M 8. To be ordered separately.</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8604 - Mounting kit, 4 countersunk bolts Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.</p> 



PRESSALIT®

Product Item number	PLUS folding seat, fixed height R5510000
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</p> <p>Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</p> <p>Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant</p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>

pressalit.com