




Product Item number	PLUS shower seat 450, manually height and sideways adjustable R7440182000
Colour Product description	<div> <div>White</div>  </div> <p>PLUS shower seat 450, with backrest. Manually height adjustable by 360 mm and sideways adjustable. The seat is foldable. Pale grey cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For fitting to horizontal PLUS wall track (not included).</p> <p>PLUS shower seat 450, with backrest. Manually height adjustable by 360 mm and sideways adjustable. The seat is foldable. Pale grey cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For fitting to horizontal PLUS wall track (not included).</p>
Adjustability Adjustment height Product capacity	Height and sideways adjustable 360 mm Maximum user mass: 200 kg. When the seat is loaded with 200 kg, each wall fastener is loaded with 23 kg. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy Features	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p> <p>Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel.</p> <p>The seat can fold up when not in use or to maximise floor space. The seat and backrest can be adjusted simultaneously thanks to a safe and easy click system.</p>

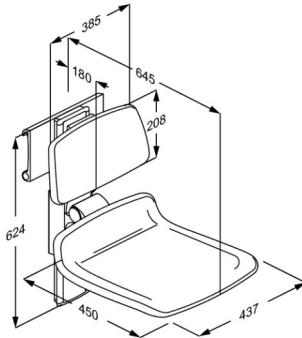

PRESSALIT®

Product	PLUS shower seat 450, manually height and sideways adjustable
Item number	R7440182000
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>PLUS shower seat.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com

PRESSALIT®