




PRESSALIT®

Product	PLUS backrest, fixed height
Item number	R6510000
Colour	<p>VVS:777538200 NRF:6190615 EAN:5708590314303 HMI: 71810 NBS:N13:426 FDA:INN/D135865</p> <p>White</p>
Product description	 <p>PLUS backrest 385 mm, soft feel. For toilet without cistern, fixed mounting. Materials: Aluminium, polyurethane foam with a plywood core. Incl. mounting kit for brick wall.</p> <p>PLUS backrest 385 mm, soft feel. For toilet without cistern, fixed mounting. Materials: Aluminium, polyurethane foam with a plywood core. Incl. mounting kit for brick wall.</p>
Adjustability	Fixed mounting
Product capacity	100 kg.
Materials	The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer.
Ergonomy	Soft and slip-resistant material that has a warm feel.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

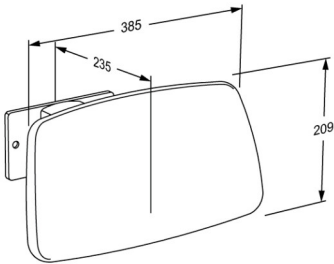

PRESSALIT®

Product	PLUS backrest, fixed height
Item number	R6510000
Cleaning	<p>PLUS backrest.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Erisan OXY+ Tvättsprit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com

PRESSALIT®