




PRESSALIT®

Product Item number	PLUS backrest, height and sideways adjustable R6530112000
Colour	VVS:777538600 NRF:6190622 EAN:5708590314488 HMI: 71808 NBS:N13:426 FDA:INN/D135865 White
Product description	 <p>PLUS backrest 385 mm, soft feel. Manually height adjustable by 150 mm and sideways adjustable. Anthracite cover plate for vertical wall track. Materials: Aluminium, polyurethane foam with a plywood core. For fitting to horizontal PLUS wall track (not included).</p> <p>PLUS backrest 385 mm, soft feel. Manually height adjustable by 150 mm and sideways adjustable. Anthracite cover plate for vertical wall track. Materials: Aluminium, polyurethane foam with a plywood core. For fitting to horizontal PLUS wall track (not included).</p>
Adjustability Adjustment height Product capacity	Height and sideways adjustable 150 mm 100 kg.
Materials Ergonomy	<p>The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer.</p> <p>Soft and slip-resistant material that has a warm feel.</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

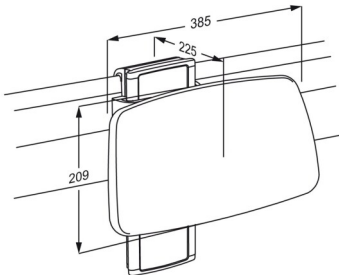

Product	PLUS backrest, height and sideways adjustable
Item number	R6530112000
Cleaning	<p>PLUS backrest.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Erisan OXY+ Tvättsprit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>