



PRESSALIT®

| | |
|--------------------------------------|---|
| Product Item number | PLUS grab bar, corner with vertical grip 30" x 17.9" x 42.9" RT138112 |
| Color | VVS:777568712 NRF:6311234 RSK:8980591 EAN:5708590318790 HMI: 25429 NBS:N13:436 Anthracite gray |
| Product description | <p>PLUS grab bars are high-quality products and they can be combined in countless ways and they provide the final touch to furnishing a functional and flexible bathroom. A reliable solution with many different options, which can easily lay the grounds for increased independence.</p> <p>PLUS grab bar combination, 1.3", reversible right/left, A:30" x B:17.9" x C: 42.9". Material: Powder-coated aluminium and plastic. Excl. mounting kit.</p> |
| Product capacity | Maximum load: 440 lbs. When the grab bar is loaded with 440 lbs, each wall fastener is loaded with 123 lbs. <i>Grab bars are tested with a 1.5 safety factor in accordance with the requirements in ISO 17966:2016.</i> |
| Materials | Core material: Aluminium Surface material: Polyester lacquer Cover cap: Cicoloy plastic Strap: Nylon (PA) |
| Ergonomy | The 1.3" diameter enables the rail to be held comfortably and firmly. The projection of 2" from the wall makes it easy to take hold of the rail. |
| Warranty | Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. |



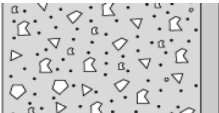



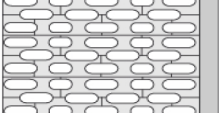

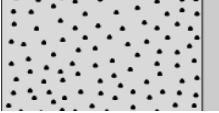
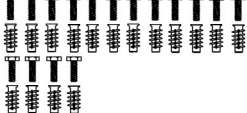
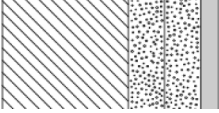
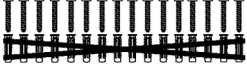
PRESSALIT®

| | |
|----------------------------------|--|
| Product | PLUS grab bar, corner with vertical grip 30" x 17.9" x 42.9" |
| Item number | RT138112 |
| Cleaning | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p> |
| Declaration of conformity | <p>This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.</p> |



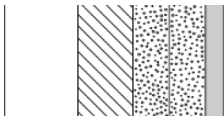
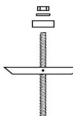

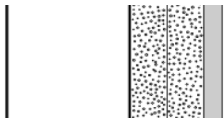
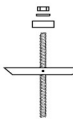

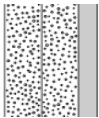

| | |
|--------------------|--|
| Product | PLUS grab bar, corner with vertical grip 30" x 17.9" x 42.9" |
| Item number | RT138112 |

Mounting in different wall types

| | |
|---|---|
|  | <p>Concrete >B25</p> <p>V8618 - Mounting kit, 16 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>  |
|  | <p>Solid Brick >= Mz12</p> <p>V8618 - Mounting kit, 16 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>  |
|  | <p>Solid Sand Lime Brick >= Hlz12</p> <p>V8618 - Mounting kit, 16 bolts Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>  |
|  | <p>Aerated Concrete P4</p> <p>V8668 - Mounting kit, 16 bolts M5 x 50 - plugs FTP K 6, for aerated concrete. To be ordered separately.</p>  |
|  | <p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8618 - Mounting kit, 16 bolts Ø5 x 50 - plugs SX8. NB! Plugs are not needed. To be ordered separately.</p>  |



PRESSALIT®

| Product Item number | PLUS grab bar, corner with vertical grip 30" x 17.9" x 42.9" RT138112 | | |
|---|---|---|---|
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (16 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (16 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. |  | |

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®