




PRESSALIT®

Product Item number	PLUS backrest, fixed height R6515277
Colour Product description	<div>  </div> Lime green PLUS backrest, 385 mm, soft feel, for SELECT toilet lifter. Materials: Aluminium, polyurethane foam with a plywood core. For mounting directly on SELECT toilet lifter, fixing kit is included. PLUS backrest, 385 mm, soft feel, for SELECT toilet lifter. Materials: Aluminium, polyurethane foam with a plywood core. For mounting directly on SELECT toilet lifter, fixing kit is included.
Adjustability Product capacity	Fixed mounting 100 kg.
Materials Ergonomy	The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer. Soft and slip-resistant material that has a warm feel.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

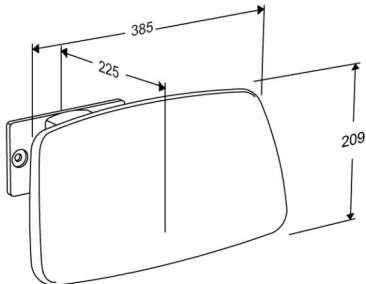

Product	PLUS backrest, fixed height
Item number	R6515277
Cleaning	<p>PLUS backrest</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Erisan OXY+</p> <p>Tvättspit 70%</p> <p>CaviCide surface disinfectant decontaminant cleaner</p> <p>Rodalon Desinfektionsspray</p> <p>Terralin protect (SCHÜLKE+)</p> <p>Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS:</p> <p>Wetrok Caledor (Wetrok)</p> <p>Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>