




Product	PLUS replacement shower seat 450 with aperture, manually height adjustable
Item number	R7401000 VVS:777534200 EAN:5708590328430
Colour	White
Product description	<p>PLUS shower seat 450 with aperture, Manually height adjustable by 230 mm. The seat is foldable.</p> <p>Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile.</p> <p>For fitting to vertical PLUS wall track (not included).</p>
Adjustability	Height adjustable
Adjustment height	230 mm
Product capacity	<p>Maximum user mass: 200 kg.</p> <p>When the seat is loaded with 200 kg, each wall fastener is loaded with 47 kg.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>
Materials	The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.
Ergonomy	Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel.
Features	The seat can fold up when not in use or to maximise floor space.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



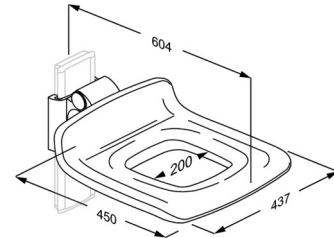

PRESSALIT®

Product	PLUS replacement shower seat 450 with aperture, manually height adjustable
Cleaning	<p>PLUS shower seat. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
Declaration of conformity	<p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007.</p>



www.pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®