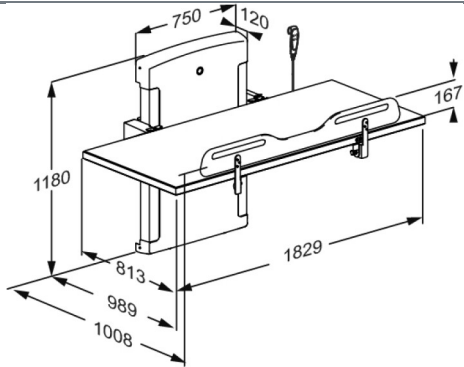




<p>Product Item number</p>	<p>CT 4100 change table, electrically foldable, electrically height adjustable R8585572000</p>
<p>Colour Product description</p>	<p>VVS: 777541200 EAN: 5708590358147 HMI: 125837</p> <p>White</p>  <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>The CT 4000 series has been designed for public changing facilities and schools.</p> <p>CT 4100 change table, electrically foldable. Electrically height adjustable 700 mm; from 300 mm -1000 mm, with buttons placed in table front. Incl. safety rail, body belt and wired remote control. Lying surface 813 x 1829 mm, in compact laminate. Outer measurements 1008 mm x 1829 mm. Excl. mounting kit.</p>
<p>Adjustability Adjustment height Product capacity Speed upwards Noise level</p>	<p>Height adjustable</p> <p>700 mm</p> <p>Maximum user weight (MUM): 227 kg. Safe working load (SWL): 232 kg. When the table is loaded with 227 kg, each mounting bolt supports 37,5 kg <i>This product has been tested statically with a load of 580 kg as the IEC 60601-1 standard requires testing with safety factor 2.5.</i></p> <p>Speed up mm: 15 mm/s Speed down mm: 20 mm/s</p> <p>48 dB</p>



Product	CT 4100 change table, electrically foldable, electrically height adjustable
Item number	R8585572000
Materials	Powder-laquered stainless steel, lying surface in compact laminate.
Ergonomy	The ergonomically designed safety rail has a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user. Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.
Mains connection (V)	230
Actuator (motor)	LA40 Careline (Linak) 4000 N, IPX6
Control box	CA63 Careline (Linak), IPX6
Wired remote control	HB30 (Linak) IP66
Power consumption (maximum load)	Power consumption by maximum load, watt: 420 W Power consumption by standby: 2 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date.

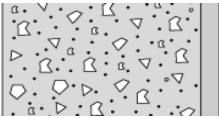
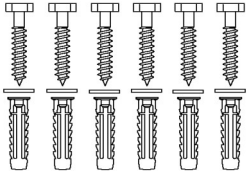

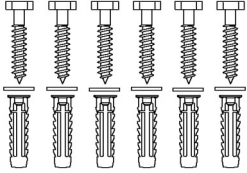
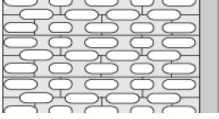
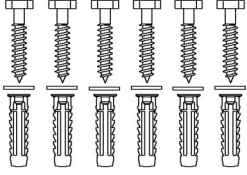
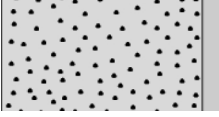
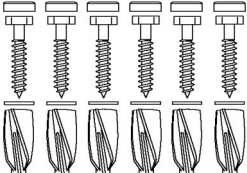
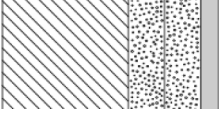
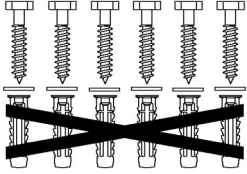


Product	CT 4100 change table, electrically foldable, electrically height adjustable
Item number	R8585572000
Cleaning	<p>CT 4000, 4100, 5000 The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningsprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
WEEE	Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.
Declaration of conformity	<p>CE UK MDR-MASK-RoHS This product is in compliance with all relevant requirements of:</p> <ul style="list-style-type: none">- Regulation (EU) 2017/745 on medical devices- Council Directive 2006/42/EC on machinery- Council Directive 2011/65/EC on RoHS. For the UK:- Medical Device Regulation (MDR) 2002 on medical devices.- Supply of Machinery (Safety) Regulations 2008- Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012.

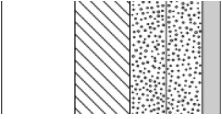
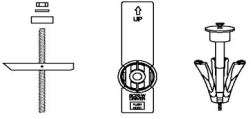
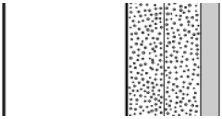
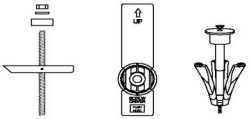
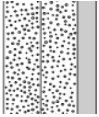



Product	CT 4100 change table, electrically foldable, electrically height adjustable
Item number	R8585572000

Installation in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.	
	Aerated Concrete P4 V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.	



Product Item number	CT 4100 change table, electrically foldable, electrically height adjustable R8585572000	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.	

It is the responsibility of the installer to ensure that suitable wall fasteners are used.