

Product	PLUS shower seat 450, manually height adjustable				
Item number	R7434035				
	NRF:6317123 EAN:5708590359786 HMI: 124335 FDA: INN/D135865				
Colour	Red (RAL 3003)				
Product description	PLUS shower seats are offered in a number of colours. They are available as both fixed height, electrically adjustable and manually adjustable, with and without aperture, and with an option of back and armrests. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.				
	PLUS shower seat 450, with back- and armrests. Manually height adjustable by mm, without the use of tools. Armrests are independently height adjustable. The and the armrests are foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluit profile. The non-slip gripping surface of the armrests are made of soft SEBS incorporating a glass-reinforced nylon core. Excl. mounting kit. Replaces R7430				
Adjustability	Height adjustable				
Adjustment height	195 mm				
Product capacity	Maximum user mass (seat): 200 kg. Maximum load for armrests is 55 kg. When the seat is loaded with 200 kg, each wall fastener is loaded with 24 kg. The seat has been tested statically with safety factor 1.5, which is a requirement				
	according to the applied assistive products standards.				



Product	PLUS shower seat 450, manually height adjustable			
Item number	R7434035			
Materials	The shower seat is made of polyurethane (PU) foam with a plywood core and this covered with an in-mould coated surface that is protected using a clear PU lacque. The backrest is made of polyurethane (PU) foam protected using a clear PU lacque. The arms are primarily constructed of 10-micrometre anodized aluminium, with a slip gripping surface at the end made of soft SEBS and inforporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal component fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.			
Ergonomy	Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an eand comfortable surface to lean on and push against.			
Features	The seat, backrest and armrests can be adjusted simultaneously thanks to a safe a easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximise floor space.			
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.			



Product	PLUS shower seat 450, manually height adjustable
Item number	R7434035
Cleaning	PLUS shower seat. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.
	The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.
	DISINFECTANTS: Erisan OXY+ Tvättsprit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)
	NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)
	Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.
	The product must be wiped down with a soft rag wrung out in clear water immediatel after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.
	DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the part must be wiped over with a damp cloth wrung out in clear water.
Declaration of conformity	CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices.



Product		
Item number		

PLUS shower seat 450, manually height adjustable

R7434035


Installation in different wall types


a. V
C
D C C C C C
D. C. O. V.
\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\
1. 4

Concrete >B25

V8612 - Mounting kit, 10 coach screws \emptyset 8 x 70, plugs SX10.

To be ordered separately.

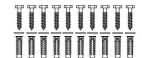


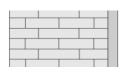


Solid Brick >= Mz12

V8612 - Mounting kit, 10 coach screws Ø8 x 70, plugs SX10.

To be ordered separately.

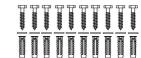


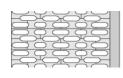

Solid Sand Lime Brick >= Hlz12

V8612 - Mounting kit, 10 coach screws \emptyset 8 x 70, plugs SX10.

To be ordered separately.

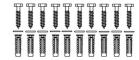





Aerated Concrete P4

V8662 - Mounting kit, 10 countersunk bolts M8 x 35, plugs FTP M 8.

To be ordered separately.

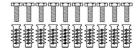


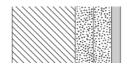

Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood

V8612 - Mounting kit, 10 coach screws Ø8 x 70, plugs \$x10

NB! Plugs are not needed. To be ordered separately.

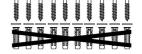



PRESSALIT®

Product Item number	PLUS shower seat 450, manually height adjustable R7434035		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	3	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	⇒ 5 3 3 3 3 3 3 3 3 3 3	

pressalit.com